**SureCC Color Compensation Kit I**

Art. No. F4009

Version 2.0

To apply a Color Compensation file for data analysis please open the multiplex assay and click the „Filter Combination“ button to select the filter combination you want to display. In the Color Comp drop-down menu, select “in Database” and choose the stored Color Compensation file you want to apply to the assay. The “Color Comp” button switches to “Color Comp (On)” to confirm that Color Compensation is applied.

Note: As the Color Compensation is instrument-specific, it is necessary to generate a Color Compensation object for every LightCycler® 480 system. A new object also has to be created after the optical system has been repaired.

**Product Information**

- Material Safety Data Sheet

**Technical Support**

For further questions please contact your distributor or send an e-mail to [info@congen.de](mailto:info@congen.de).

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In the next step please select the "Sample Editor" button on the left side. Here please choose in the Step 1 „Select Workflow“-box the "Color Comp" tab. Now please select in Step 2 dialog box the subset of the programming step before. In the "Dominant Channel" dialog fields select the corresponding dyes (FAM, VIC, ROX, Cy5) in the different wells. For the replicate samples containing no dye, select "Water". It is not necessary to give a sample name for the wells.

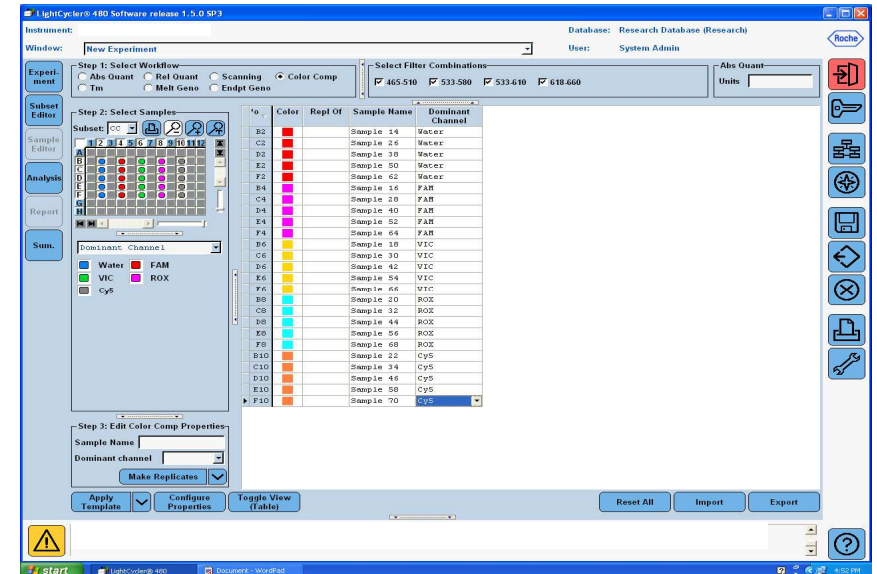
Fig. 8

Please save your changes and the experiment by clicking on the disc icon. Place the PCR plate (which contains the reactions prepared before) into the LightCycler® 480I or II. Click the "Start run" button in the Experiment window (on the left side) to start the experiment.

**Evaluation and creation of a CC-File**

After the experiment ends, click the "Analysis" button in the left side of the window to open the "Create New Analysis" dialog box. Select "Color Compensation" and confirm your choice.

In the next dialog-box please select the subset and confirm your selection. Click now the button „Calculate“ in the analysis window to perform color compensation analysis. Click „Save CC Object“ button to save the Color Compensation file in the "CCC" folder. The Color Compensation file is now ready to use.

The final screen should look like the screenshot below:

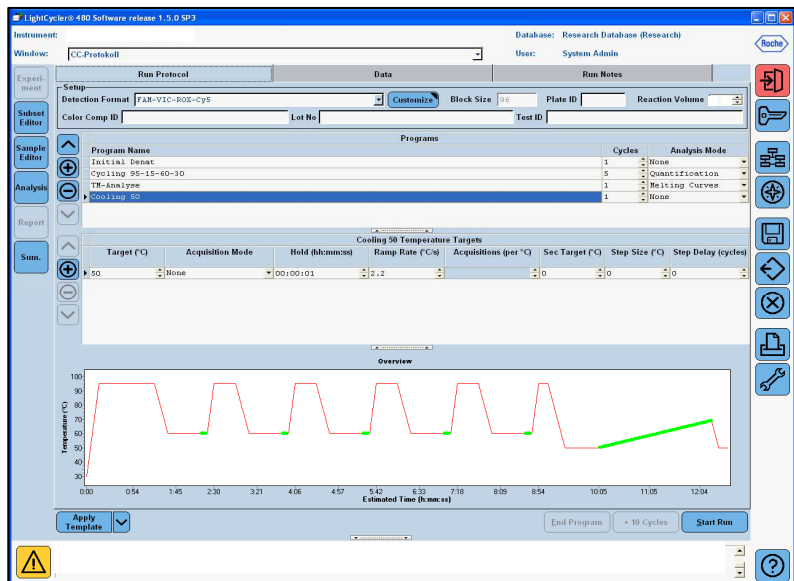
Fig. 6

For selecting the right PCR layout please click on the "Subset Editor" button on the left side. Please create a subset through a click on the Plus-Button and give the subset a name for example "CC Layout". Now please press the Strg-key and the left mouse button simultaneously. Mark all wells with the solution inside the plate. Then click on the "Apply" button and the subset is ready. The final screen should look like the screenshot below:

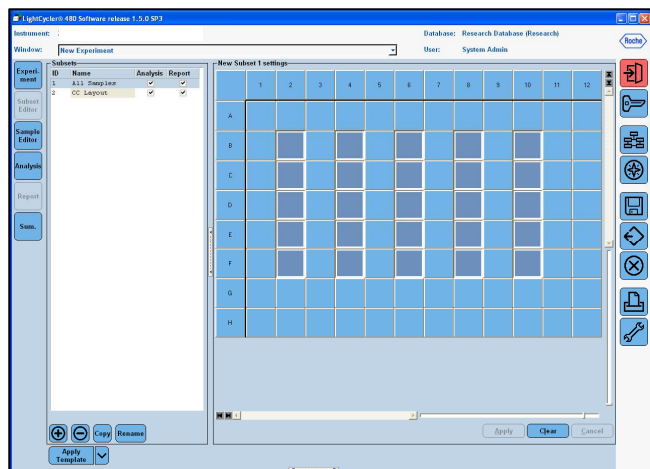
Fig 7

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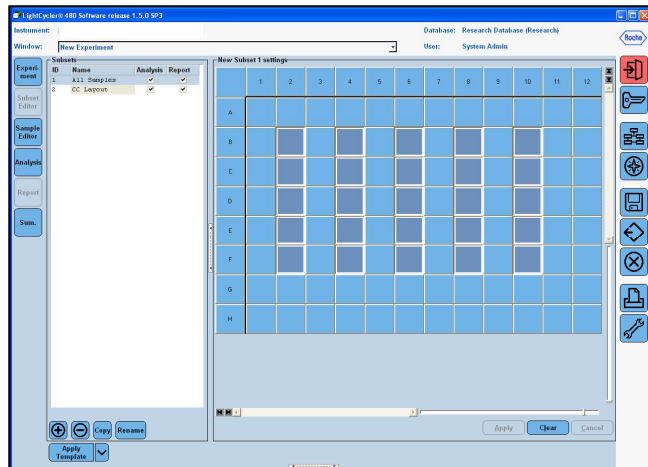
Bild 7

After the programming of the detection format please choose the button „New Experiment“ in the screen of figure 2. Following window will open:

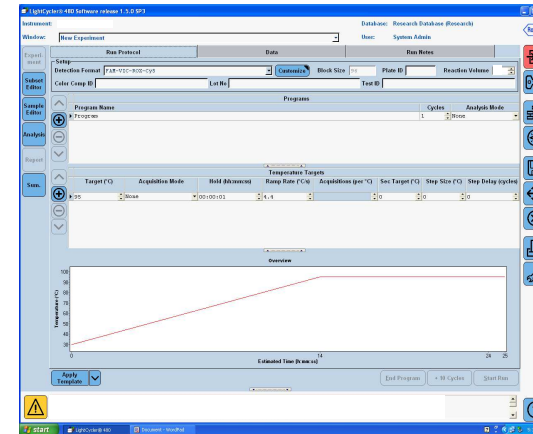
Fig. 4

As the first step select the right detection format (please choose the right name – see procedure in figures 2 and 3) and correct reaction volume (20 µl).

Following protocol steps are to program:  
(Please pay attention to the correct Cycle number und Analysis Mode.)

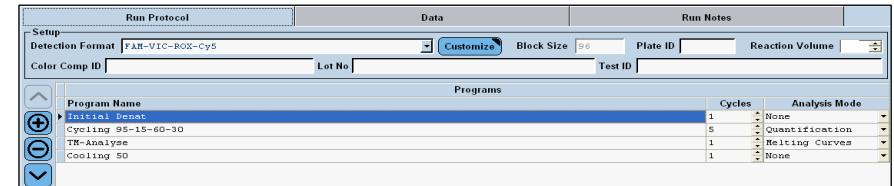
Fig. 5

Please program following protocol steps, temperature targets and measurement points:

Program	Cycles / Analysis Mode	Temperature targets			
		Target [°C]	Acquisition Mode	Hold [hh:mm:ss]	Ramp rate [°C/s]
Initial Denat.	1 / none	95	none	00:00:30	4.4
Cycling 95-15-60-30	5 / Quantification	95	none	00:00:15	4.4
		60	single	00:00:30	2.2
TM-Analyse	1 / Melting Curves	95	none	00:00:01	4.4
		50	none	00:00:30	2.2
		70	continuous		0.14 (Acquisitions per °C = 1)
Cooling 50	1 / none	50	none	00:00:01	2.2

2. Setup and start of the experiment

Open the LightCycler® 480I or II software. Please program under the button for tools the necessary detection format:

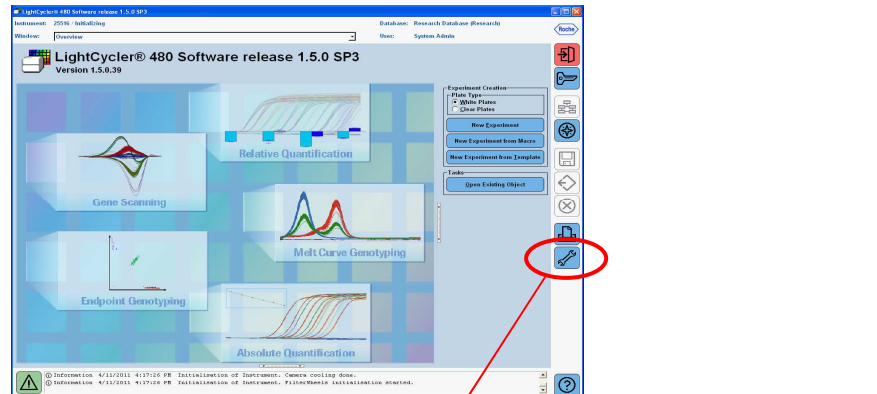
Fig. 2

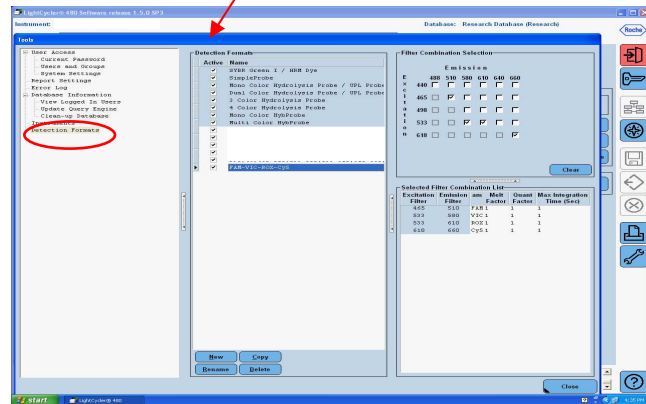
Fig. 3

Generate and save the following detection format for example as „FAM-VIC-ROX-Cy5“:

Dye-Name	Filter Combination LightCycler® 480I	Filter Combination LightCycler® 480II
FAM	450 – 500	465 - 510
VIC	523 – 568	533 - 580
ROX	558 – 610	533 - 610
Cy5	615 – 670	618 - 660

- Value for Quant- and Melt Factors as well as the Integration Time please set 1

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**Description**

This kit is intended for the generation of color compensation objects for the LightCycler® 480I and II system. With 5 solutions it is possible to perform a color compensation experiment. Fluorescence data are collected and used to generate a color compensation file containing information for correcting crosstalk between detection channels in tetraplex real-time PCR runs. The generated color compensation object can subsequently be used to analyze tetraplex real-time PCR experiments with the dyes FAM, VIC, ROX and Cy5.

**Kit components and storage**

The kit is designed for 3 calibration runs.

- 1x Blank (400 µl) **(Code 1)**
- 1x Dye 1 (400 µl) **(Code 2)**
- 1x Dye 2 (400 µl) **(Code 3)**
- 1x Dye 3 (400 µl) **(Code 4)**
- 1x Dye 4 (400 µl) **(Code 5)**

Store all reagents at -20°C and protected from light.

**Additionally required equipment and materials**

- Real-time PCR device: LightCycler® 480I or II
- Real-time PCR consumables (plates, foils)
- Pipettes with filter tips

**Protocol**

1. Preparation of the color compensation plate

For a color compensation experiment it is necessary to pipette 5 reactions with 20 µl for each dye and 5 reactions with 20 µl for the fluorescence background. See the layout below:

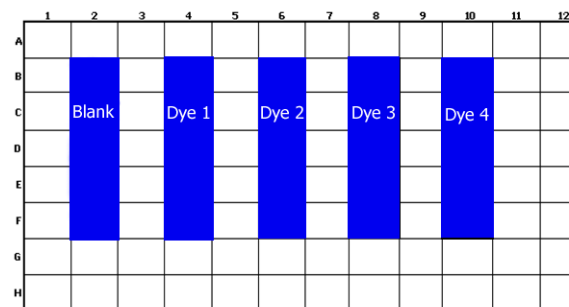
Fig. 1

Allow the reagents to thaw, mix by vortexing and centrifuge before opening and use.

Solution	Dye	Amount per reaction	Pipetting of 20 µl in following wells
Blank	Background	20.0 µl	B2, C2, D2, E2, F2
Dye 1	FAM	20.0 µl	B4, C4, D4, E4, F4
Dye 2	VIC	20.0 µl	B6, C6, D6, E6, F6
Dye 3	ROX	20.0 µl	B8, C8, D8, E8, F8
Dye 4	Cy5	20.0 µl	B10, C10, D10, E10, F10

Close the plate with an optical sealing foil after pipetting the solutions.