

AO ZON PREP[®]

Product Code: P112 / P112B

Immunoaffinity columns for use in conjunction with HPLC or LC-MS/MS.
For in vitro use only.

P112/V11/17.09.21

www.r-biopharm.com



R-BIOPHARM
RHÔNE LTD

Contents

	Page
Test Principle.....	4
Reagents Not Provided.....	4
Accessory Products.....	4
Hazards.....	4
Recommended Methods and Application Notes.....	4
Decontamination.....	5
Storage & Shelf Life.....	5
Sampling.....	5
Sensitivity.....	5
Recoveries.....	5
Column Preparation.....	5
Elution.....	6
HPLC Information.....	7
• Sample Preparation - Cereal and Baby Food.....	7
• Sample Preparation - Rye.....	8
• Preparation of Standards.....	9
• Preparation of 1,000 ng/ml Total Aflatoxin Stock Solution.....	9
• Preparation of 1,000 ng/ml Ochratoxin A Stock Solution.....	9
• Preparation of 1,000 ng/ml Zearalenone Stock Solution.....	9
• Calibration Curve.....	10
• Total Aflatoxin Standard Dilution.....	10
• Ochratoxin A Standard Dilution.....	10
• Zearalenone Standard Preparation.....	10
• Combined Working Standard.....	10
• Recommended HPLC Conditions.....	11
• Total Aflatoxin and Ochratoxin A.....	11
• Zearalenone.....	11
• HPLC Information - Typical Chromatogram.....	12
• Example HPLC Chromatogram for Maizel.....	12
• Example HPLC Chromatogram for Rye.....	13
LC-MS/MS Information.....	14
• Sample Preparation - Cereal.....	14
• Sample Preparation - Rye.....	15
• Sample Preparation - Baby Food.....	16
• Preparation of Standards.....	17
• Preparation of 1,000 ng/ml Total Aflatoxin Stock Solution.....	17
• Preparation of 1,000 ng/ml Ochratoxin A Stock Solution.....	17
• Preparation of 1,000 ng/ml Zearalenone Stock Solution.....	17
• Calibration Curve.....	18
• Total Aflatoxin Standard Dilution.....	18
• Ochratoxin A Standard Dilution.....	18
• Zearalenone Standard Preparation.....	18
• Combined Working Standard.....	18
• Recommended LC-MS/MS Conditions.....	19
• LC-MS/MS Information - Typical Chromatogram.....	20
• Example LC-MS/MS Chromatogram for Baby Food.....	20
• Example LC-MS/MS Chromatogram for Rye.....	21
Quality.....	22
Technical Support.....	22
Acknowledgement.....	22
Warranty.....	22

Test Principle

The procedure is based on monoclonal antibody technology, which makes the test highly specific, sensitive, rapid and simple to perform.

The columns contain a gel suspension of monoclonal antibody specific to the toxins of interest. Following extraction of the toxins the sample extract is filtered, diluted and passed slowly through the immunoaffinity column. Any toxins which are present in the sample are retained by the antibody within the gel suspension. The column is washed to remove unbound material and the toxins are then released from the column following elution with solvent. The eluate is collected prior to analysis by HPLC or LC-MS/MS. Aflatoxins are required to be derivatised when analysed by HPLC.

The total extraction and clean-up time takes approximately 20 minutes to perform. The result is improved clean-up and concentration of the toxins from food and feed samples giving a much cleaner chromatogram and therefore providing more accurate and sensitive detection. The columns also have the added advantage that they can be automated for large scale analysis of samples.

Reagents Not Provided

- Distilled / Deionised Water (suitable for use with HPLC, e.g. MilliQ)
- Solvents (HPLC Grade Methanol and Acetonitrile)
- Phosphate Buffered Saline (PBS) (RP202)*
- Mycotoxin Standards (Please refer to Preparation of Standards section)
- Sodium Chloride
- Sodium Hydroxide (to pH filtrate if required)
- Nitric Acid (only required when derivatising with a KOBRA® CELL)
- Potassium Bromide (only required when derivatising with a KOBRA® CELL)

Accessory Products

- Whatman No. 113 or No. 4 Filter Paper
- KOBRA® CELL (K01)*
- Immunoaffinity Column Rack (CR1)*
- Immunoaffinity Column Accessory Pack (AP01)*

* Available from R-Biopharm. Please contact your local R-Biopharm distributor for further information.

Hazards

Mycotoxins are very hazardous substances. Only laboratories equipped to handle toxic materials and solvents should perform analyses. Suitable protective clothing, including gloves, safety glasses and lab coats should be worn throughout the analysis.

Flammable solvents should be stored in an explosion-proof cabinet. Use a chemical hood and protective equipment as applicable.

Contact your local R-Biopharm distributor for a Material Safety Data Sheet for further information if required.

Recommended Methods and Application Notes

Methods are available for all matrices covered by legislation as well as additional commodities. Deviation from the methods described in our Instructions For Use and Application Notes may not achieve optimum results. Please contact your local R-Biopharm distributor for further information.

Decontamination

Prior to disposal, excess standard solutions should be treated with at least one-tenth their volume of 5 % sodium hypochlorite. Labware and contaminated waste should be immersed in 5 % sodium hypochlorite solution for 30 minutes followed by the addition of 5 % acetone for 30 minutes. Flush with copious amounts of water before disposal. After decontamination labware should be thoroughly washed. Incinerate waste if regulations permit.

Storage & Shelf Life

The columns expire 18 months from date of manufacture if stored at 2 - 8 °C or 12 months from date of manufacture if stored at 21 - 25 °C. Do not freeze.

Ensure the column has not dried out and contains buffer above the gel. It is important to note the antibody included in the immunoaffinity column can be denatured by extreme temperature or pH change.

Sampling

A representative sample should be obtained by following one of the officially recognised sampling procedures. It is recommended that a minimum of 1 kg of representative sample is finely ground and a portion (5 - 50 g dependent on method used) of this is removed and extracted.

Sensitivity

The sensitivity is dependent on the final detection system employed by the analyst. However the test sensitivity may be improved if required by increasing the volume of sample passed through the immunoaffinity column. Please note the ratio of solvent to phosphate buffered saline (PBS) should be maintained.

Recoveries

If an analyst wishes to account for losses during extraction it is recommended a spiked sample of the same commodity type as the material being tested is analysed following the complete procedure as a reference standard. The recoveries obtained with the spiked sample can be used to correct the results obtained with the test sample.

Column Preparation

Immunoaffinity columns should be at ambient temperature before use. Remove the cap from the top of the column and discard. Firmly attach the column to a glass syringe barrel using an adapter and place in an immunoaffinity column rack or clamp stand.

Elution


In order to fully elute the toxin/s from the immunoaffinity column it is vital that the solvent is in contact with the antibody within the gel suspension for a sufficient period of time. This ensures that all of the bonds between the antibody and the toxin are broken, ultimately releasing all of the toxin from the column for analysis with the detection system of choice

To ensure that the solvent is in contact with the antibody gel for a sufficient period of time any of the following elution methods can be used: -

Backflushing (this is the preferred method of choice at R-Biopharm): backflush by gently raising and lowering the syringe plunger during passage of the solvent through the column. This process will reverse the direction of flow of the eluate through the gel. This should be repeated 3 times before collecting the eluate. Proceed to the next step in the method.

Application of small volumes of solvent: apply the volume of solvent required for elution in two or three smaller aliquots. Allow each aliquot to remain in contact with the gel suspension for a minimum of 30 seconds before allowing each to pass fully through the gel suspension for collection. Proceed to the next step in the method.

Incubation with solvent: apply the full volume of solvent required for elution and allow 2-3 drops of the solvent to pass through the column for collection. Allow the remainder of the solvent to remain in contact with the gel suspension for a minimum of 60 seconds before allowing it to pass through the gel suspension for collection. Proceed to the next step in the method.

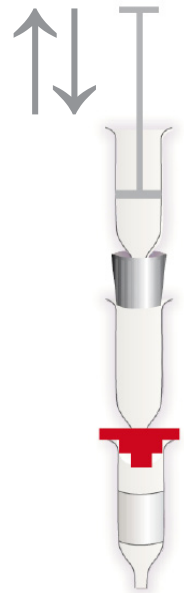

HPLC Information

• Sample Preparation - Cereal and Baby Food

This method has been tested on a number of commodities including wheat, maize and baby food.

1. Weigh 25 g of ground sample and 5 g of sodium chloride into a 1 litre capacity, solvent resistant blender jar.
2. Add 100 ml of 80 % methanol and blend at high speed for 2 minutes.
3. Filter the sample through Whatman No. 113 or No. 4 filter paper, or centrifuge at 4,000 rpm for 10 minutes.
4. Dilute 10 ml of the filtrate with 40 ml of phosphate buffered saline (PBS).
5. Pass 20 ml of the filtrate (equivalent to 1 g of sample) through the column at a flow rate of 2 ml per minute (or the sample can be allowed to pass through the column by gravity if preferred). A slow, steady flow rate is essential for the capture of the toxins by the antibody.
6. Wash the column by passing 20 ml of PBS through at a flow rate of approximately 5 ml per minute. Pass air through the column to remove residual liquid.
7. Elute the toxins from the column at a flow rate of 1 drop per second using 1.5 ml of 100 % methanol and collect in an amber glass vial. Please refer to the Elution section for further information.
8. Following elution pass 1.5 ml of water through the column and collect in the same vial to give a 3 ml total volume.
9. Pipette 1 ml of eluate into two separate amber glass vials and inject 100 μ l onto each HPLC system for analysis of zearalenone and combined analysis of total aflatoxin and ochratoxin A.

HPLC Information

• Sample Preparation - Rye

1. Weigh 25 g of ground sample into a 1 litre capacity, solvent resistant blender jar.
2. Add 100 ml of 60 % acetonitrile and blend at high speed for 2 minutes.
3. Filter the sample through Whatman No. 113 or No. 4 filter paper, or centrifuge at 4,000 rpm for 10 minutes.
4. Dilute 10 ml of the filtrate with 40 ml of phosphate buffered saline (PBS).
5. Filter the diluted extract through glass microfibre filter paper.
6. Pass 20 ml of the filtrate (equivalent to 1 g of sample) through the column at a flow rate of 2 ml per minute (or the sample can be allowed to pass through the column by gravity if preferred). A slow, steady flow rate is essential for the capture of the toxin by the antibody.
7. Wash the column by passing 20 ml of PBS through at a flow rate of approximately 5 ml per minute. Pass air through the column to remove residual liquid.
8. Elute the toxins from the column at a flow rate of 1 drop per second using 1.5 ml of 100 % methanol and collect in a 5 ml amber glass vial. Please refer to the Elution section for further information.
9. Following elution pass 1.5 ml of water through the column and collect in the same vial to give a 3 ml total volume.
10. Pipette 1 ml of eluate into two separate amber glass vials and inject 100 μ l onto each HPLC system for analysis of zearalenone and combined analysis of total aflatoxin and ochratoxin A.

HPLC Information

Preparation of Standards

- **Aflatoxin Stock Solution**

It is advised to start with a 1,000 ng/ml total aflatoxin stock solution.

Note: The ratio of B1, B2, G1 and G2 may vary in each standard. Please note the correct ratio for the standard purchased.

- **Ochratoxin A Stock Solution**

It is advised to start with a 1,000 ng/ml ochratoxin A stock solution.

- **Zearalenone Stock Solutions**

It is advised to start with a 1,000 ng/ml zearalenone stock solution.

HPLC Information

Calibration Curve

• Total Aflatoxin Standard Dilution

1. Measure 1 ml of 100 % methanol into an amber vial.
2. Remove 100 µl to waste.
3. Add 100 µl of 1,000 ng/ml total aflatoxin standard to give a 100 ng/ml total aflatoxin solution.

• Ochratoxin A Standard Dilution

1. Measure 1 ml of 100 % methanol into an amber vial.
2. Remove 100 µl to waste.
3. Add 100 µl of 1,000 ng/ml ochratoxin A standard to give a 100 ng/ml ochratoxin A solution.

• Zearalenone Standard Dilution

1. 1,000 ng/ml zearalenone standard solution should be used, no dilution is required.

• Combined Working Standard

It is recommended to run at least a 3 - 6 point calibration curve. In constructing a suitable curve the levels of the calibration standards should bracket or include the range of expected results. The diluted standard solutions should be prepared fresh on the day of analysis and used within a 24 hour period.

Example of how to prepare a four point calibration curve (can be modified according to legislative requirements or contamination levels):

1. Standard 4:
 - Measure 3 ml of 100 % methanol into an amber vial and discard 780 µl to waste.
 - Add 240 µl of 100 ng/ml total aflatoxin standard, 240 µl of 100 ng/ml ochratoxin A standard and 300 µl of 1,000 ng/ml of zearalenone standard.
 - Add 3 ml of water (equivalent to 4 ng/ml total aflatoxin, 4 ng/ml ochratoxin A and 50 ng/ml zearalenone). Vortex for 20 seconds.
2. Standard 3: Take 3 ml of standard 4 and add 3 ml of 50 % methanol (equivalent to 2 ng/ml total aflatoxin, 2 ng/ml ochratoxin A and 25 ng/ml zearalenone).
3. Standard 2: Take 3 ml of standard 3 and add 3 ml of 50 % methanol (equivalent to 1 ng/ml total aflatoxin, 1 ng/ml ochratoxin A and 12.5 ng/ml zearalenone).
4. Standard 1: Take 3 ml of standard 2 and add 3 ml of 50 % methanol (equivalent to 0.5 ng/ml total aflatoxin, 0.5 ng/ml ochratoxin A and 6.25 ng/ml zearalenone).
5. Dispense 1 ml of each standard into 2 amber vials and inject 100 µl of each solution onto the HPLC system.

HPLC Information

Recommended HPLC Conditions

- Total Aflatoxin and Ochratoxin A

HPLC Conditions

Derivatisation	KOBRA® CELL at 100 µA setting		
Guard Cartridge	Inertsil ODS-3 5 µm, 4 mm x 10 mm or (Hichrom) equivalent		
Analytical Column	Inertsil ODS-3V 5 µm, 4.6 mm x 150 mm (Hichrom) or equivalent		
Mobile Phase	Solution A: Water : Methanol (55 : 45 v/v) Solution B: Water : Methanol (20 : 80 v/v) Add 119 mg of potassium bromide and 350 µl 4 M Nitric Acid to 1 litre of mobile phase A and B. Prepare fresh on day of analysis.		
Gradient Conditions	Time (min)	% Solution A	% Solution B
	0	100	0
	14	100	0
	16	35	65
	30	35	65
	31	100	0
	40	100	0
HPLC Pump	From preferred supplier		
Flow Rate	0.8 ml per minute		
Fluorescence Detector	Time (min)	Excitation (nm)	Emission (nm)
	0	365	442
	17	333	463
Column Heater	Maintain guard and analytical columns at 40 °C		
Integrator / Data Control System	From preferred supplier		
Injector	Autosampler / Rheodyne valve		
Injection Volume	100 µl		
Elution Order	G2, G1, B2, B1, ochratoxin A		



- Zearalenone

HPLC Conditions

Guard Cartridge	Inertsil ODS-3 5 µm, 4 mm x 10 mm or (Hichrom) equivalent		
Analytical Column	Inertsil ODS-3V 5 µm, 4.6 mm x 150 mm (Hichrom) or equivalent		
Mobile Phase	Acetonitrile : Water : Methanol (46 : 46 : 8 v/v) Prepare fresh on day of analysis.		
HPLC Pump	To deliver mobile phase		
Flow Rate	1.0 ml per minute		
Fluorescence Detector	Excitation: 274 nm		
	Emission: 455 nm		
Column Heater	Maintain guard and analytical columns at 40 °C		
Integrator / Data Control System	From preferred supplier		
Injector	Autosampler / Rheodyne valve		
Injection Volume	100 µl		



HPLC Information - Typical Chromatograms

Example HPLC Chromatogram for Maize

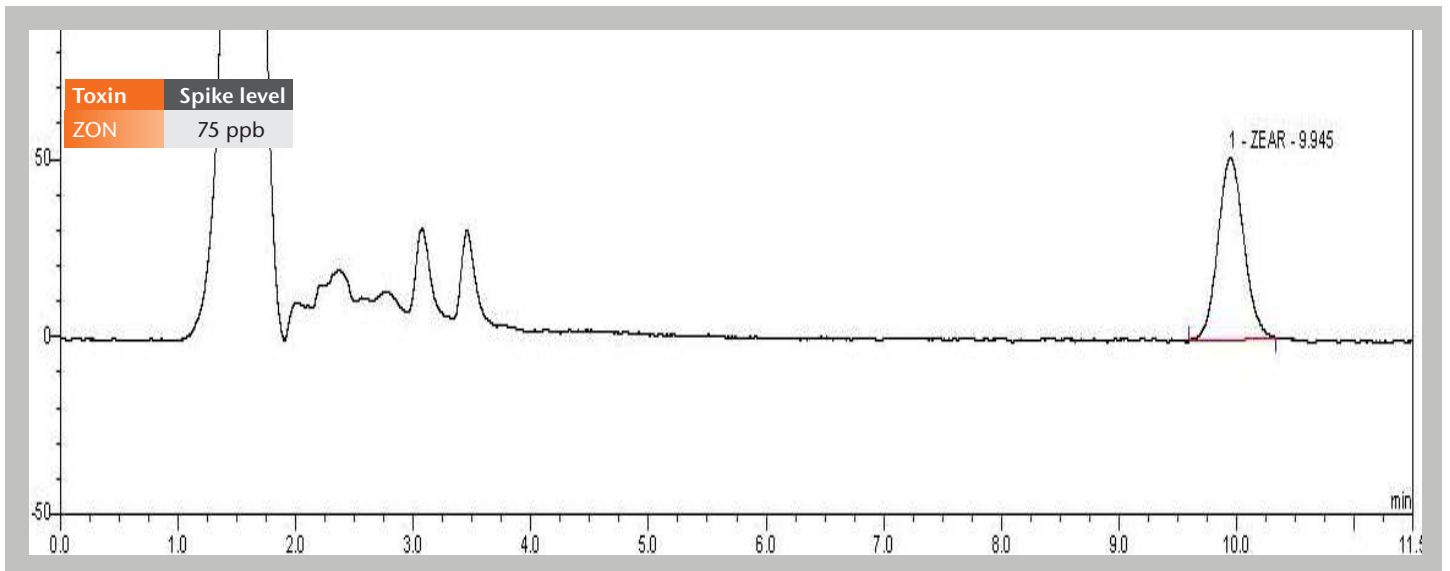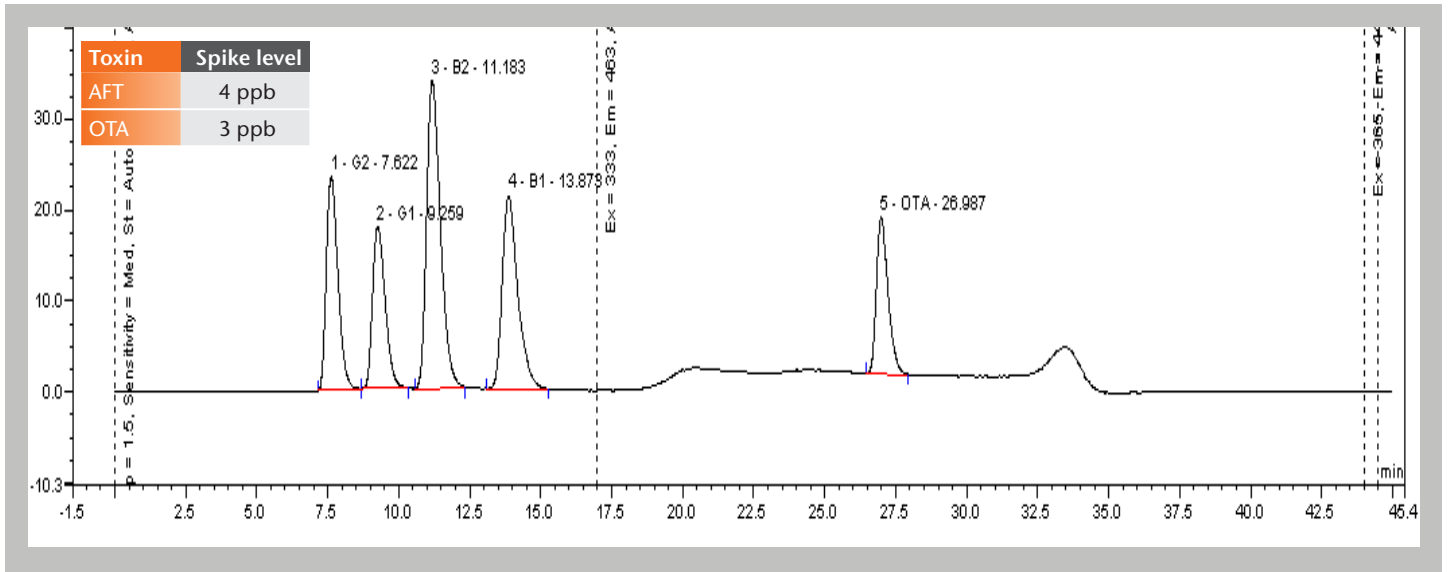

HPLC Information - Typical Chromatograms

Example HPLC Chromatogram for Rye

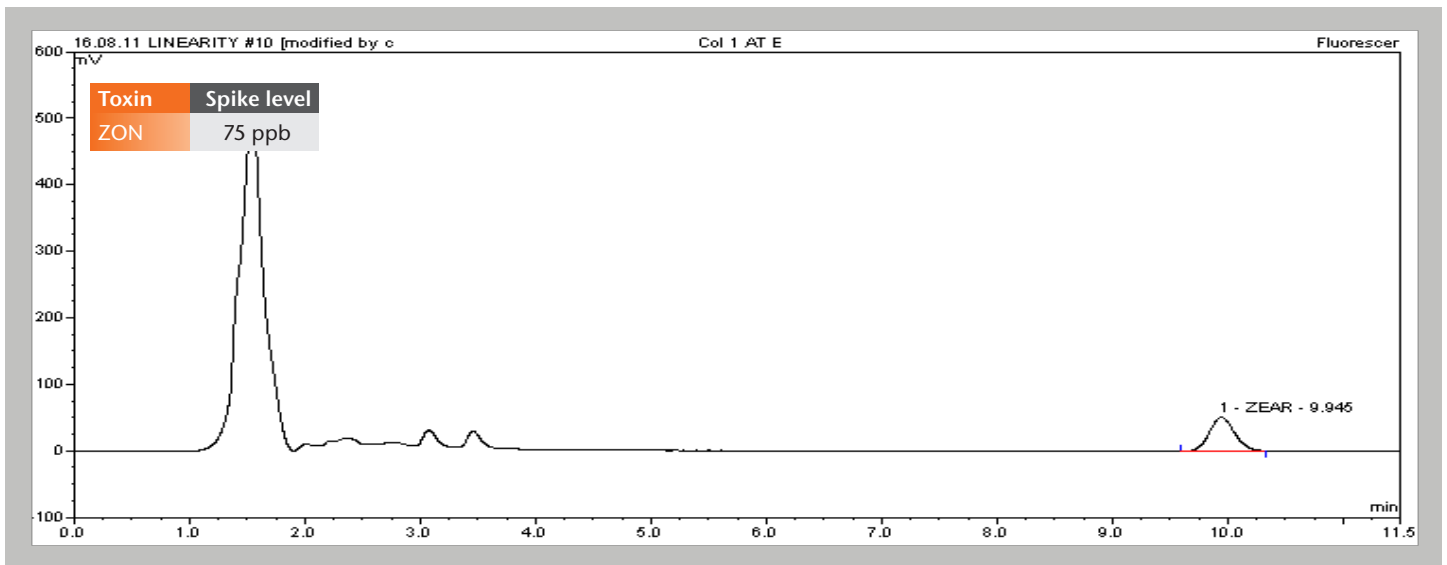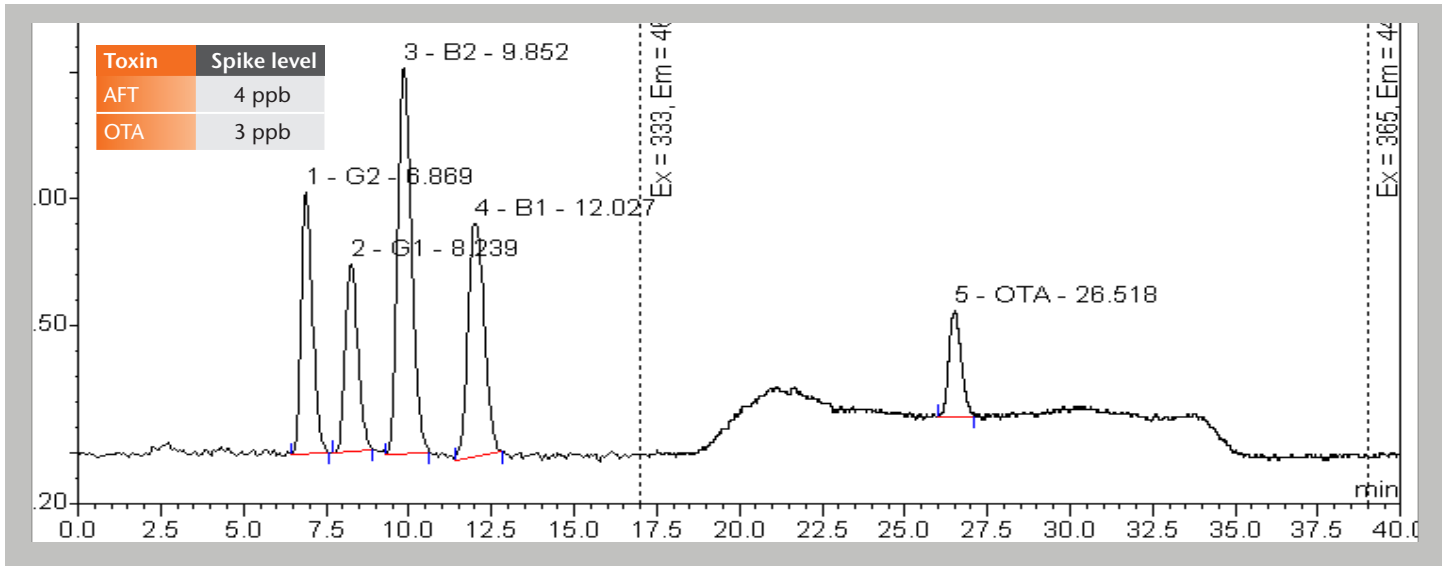

LC-MS/MS Information

• Sample preparation - Cereal

1. Weigh 25 g of ground sample and 5 g of sodium chloride into a 1 litre capacity, solvent resistant blender jar.
2. Add 100 ml of 80 % methanol and blend at high speed for 2 minutes.
3. Centrifuge at 4,000 rpm for 10 minutes.
4. Dilute 10 ml of the filtrate with 40 ml of phosphate buffered saline (PBS).
5. Filter the diluted extract through glass microfibre filter paper.
6. Pass 20 ml of the filtrate (equivalent to 1 g of sample) through the column at a flow rate of 2 ml per minute (or the sample can be allowed to pass through the column by gravity if preferred). A slow, steady flow rate is essential for the capture of the toxins by the antibody.
7. Wash the column by passing 20 ml of water through at a flow rate of approximately 5 ml per minute. Pass air through the column to remove residual liquid.
8. Elute the toxins from the column at a flow rate of 1 drop per second using 1.5 ml of 100 % methanol and collect in an amber glass vial. Please refer to the Elution section for further information.
9. Following elution pass 1.5 ml of water through the column and collect in the same vial to give a 3 ml total volume.
10. Inject 25 μ l onto the LC-MS/MS system

LC-MS/MS information

• Sample preparation - Rye

1. Weigh 25 g of ground sample into a 1 litre capacity, solvent resistant blender jar.
2. Add 100 ml of 60 % acetonitrile and blend at high speed for 2 minutes.
3. Filter the sample through Whatman No. 113 or No. 4 filter paper, or centrifuge at 4,000 rpm for 10 minutes.
4. Dilute 20 ml of the filtrate with 80 ml of phosphate buffered saline (PBS).
5. Filter the diluted extract through glass microfibre filter paper.
6. Pass 20 ml of the filtrate (equivalent to 1 g of sample) through the column at a flow rate of 2 ml per minute (or the sample can be allowed to pass through the column by gravity if preferred). A slow, steady flow rate is essential for the capture of the toxin by the antibody.
7. Wash the column by passing 20 ml of water through at a flow rate of approximately 5 ml per minute. Pass air through the column to remove residual liquid.
8. Elute the toxins from the column at a flow rate of 1 drop per second using 1.5 ml of 100 % methanol and collect in a 5 ml amber glass vial. Please refer to the Elution section for further information.
9. Following elution pass 1.5 ml of water through the column and collect in the same vial to give a 3 ml total volume.
10. Inject 25 μ l onto the LC-MS/MS system.

LC-MS/MS information

• Sample preparation - Baby food

1. Weigh 25 g of ground sample and 5 g of sodium chloride into a 1 litre capacity, solvent resistant blender jar.
2. Add 100 ml of 80 % methanol and blend at high speed for 2 minutes.
3. Centrifuge at 4,000 rpm for 10 minutes.
4. Dilute 15 ml of the filtrate with 60 ml of phosphate buffered saline (PBS).
5. Filter the diluted extract through glass microfibre filter paper.
6. Pass 45 ml of the filtrate (equivalent to 2.25 g of sample) through the column at a flow rate of 2 ml per minute (or the sample can be allowed to pass through the column by gravity if preferred). A slow, steady flow rate is essential for the capture of the toxin by the antibody.
7. Wash the column by passing 20 ml of water through at a flow rate of approximately 5 ml per minute. Pass air through the column to remove residual liquid.
8. Elute the toxins from the column at a flow rate of 1 drop per second using 1 ml of 100 % methanol and collect in a 5 ml amber glass vial. Please refer to the Elution section for further information.
9. Following elution pass 1 ml of water through the column and collect in the same vial to give a 2 ml total volume.
10. Inject 25 μ l onto the LC-MS/MS system.

LC-MS/MS Information

Preparation of Standards

- **Aflatoxin Stock Solution**

It is advised to start with a 1,000 ng/ml total aflatoxin stock solution.

Note: The ratio of B1, B2, G1 and G2 may vary in each standard. Please note the correct ratio for the standard purchased.

- **Ochratoxin A Stock Solution**

It is advised to start with a 1,000 ng/ml ochratoxin A stock solution.

- **Zearalenone Stock Solutions**

It is advised to start with a 1,000 ng/ml zearalenone stock solution.

LC-MS/MS Information

Calibration Curve

• Total Aflatoxin Standard Dilution

1. Measure 1 ml of 100 % methanol into an amber vial.
2. Remove 100 µl to waste.
3. Add 100 µl of 1,000 ng/ml total aflatoxin standard to give a 100 ng/ml total aflatoxin solution.

• Ochratoxin A Standard Dilution

1. Measure 1 ml of 100 % methanol into an amber vial.
2. Remove 100 µl to waste.
3. Add 100 µl of 1,000 ng/ml ochratoxin A standard to give a 100 ng/ml ochratoxin A solution.

• Zearalenone Standard Dilution

1. 1,000 ng/ml zearalenone standard solution should be used, no dilution is required.

• Combined Working Standard

It is recommended to run at least a 3 - 6 point calibration curve. In constructing a suitable curve the levels of the calibration standards should bracket or include the range of expected results. The diluted standard solutions should be prepared fresh on the day of analysis and used within a 24 hour period.

Example of how to prepare a four point calibration curve (can be modified according to legislative requirements or contamination levels):

1. Standard 4:
 - Measure 3 ml of 100 % methanol into an amber vial and discard 780 µl to waste.
 - Add 240 µl of 100 ng/ml total aflatoxin standard, 240 µl of 100 ng/ml ochratoxin A standard and 300 µl of 1,000 ng/ml of zearalenone standard.
 - Add 3 ml of water (equivalent to 4 ng/ml total aflatoxin, 4 ng/ml ochratoxin A and 50 ng/ml zearalenone). Vortex for 20 seconds.
2. Standard 3: Take 3 ml of standard 4 and add 3 ml of 50 % methanol (equivalent to 2 ng/ml total aflatoxin, 2 ng/ml ochratoxin A and 25 ng/ml zearalenone).
3. Standard 2: Take 3 ml of standard 3 and add 3 ml of 50 % methanol (equivalent to 1 ng/ml total aflatoxin, 1 ng/ml ochratoxin A and 12.5 ng/ml zearalenone).
4. Standard 1: Take 3 ml of standard 2 and add 3 ml of 50 % methanol (equivalent to 0.5 ng/ml total aflatoxin, 0.5 ng/ml ochratoxin A and 6.25 ng/ml zearalenone).
5. Dispense 1 ml of each standard into an amber vial and inject 25 µl of each solution onto the LC-MS/MS system.

LC-MS/MS Information

• Recommended LC-MS/MS Conditions

LC Conditions

Analytical Column	Phenomenex Gemini C18, 110 Å, 150 x 3 mm, 5 µm		
Mobile Phase	Solution A: 1 mM Ammonium Formate + 0.1 % Formic Acid in Water : Methanol (95 : 05 v/v) Solution B: 1 mM Ammonium Formate + 0.1 % Formic Acid in Water : Methanol (02 : 98 v/v)		
Gradient Conditions	Time (min)	% Solution A	% Solution B
	0	80	20
	0.1	80	20
	10	10	90
	15	10	90
	15.1	80	20
	20	80	20
HPLC Pump	To deliver mobile phase		
Flow Rate	0.3 ml per minute		
Column Heater	Maintain guard and analytical columns at 40 °C		
Integrator / Data Control System	From preferred supplier		
Injection Volume	25 µl		

Mass Spectrometry Conditions


Instrument	Sciex QTRAP 4500
Mode	Multiple Reaction Monitoring (MRM)
Source Temperature	450 °C
Ion Spray	3500 V
IonSource Gas 1	50 psi
IonSource Gas 2	50 psi
Curtain Gas	50 psi

Instrument Setting

Toxin	Time Segment (min)	Precursor Ion (m/z)	Product Ions (m/z)	Dwell Time (s)	Collision Energy (V)	Cell Exit Potential (V)
AFT G2	8.81	331.1 [M+H] ⁺	256.90 (Quantifier) 189.10 (Qualifier)	20	40.90 54.40	20.00 14.00
AFT G1	8.81	329.1 [M+H] ⁺	243.00 (Quantifier) 199.90 (Qualifier)	20	36.60 53.70	15.00 15.00
AFT B2	9.24	315.2 [M+H] ⁺	287.10 (Quantifier) 259.00 (Qualifier)	20	34.90 38.60	17.00 17.00
AFT B1	9.59	313.2 [M+H] ⁺	285.10 (Quantifier) 241.00 (Qualifier)	20	30.70 48.20	17.00 17.00
OTA	12.74	404.1 [M+H] ⁺	239.00 (Quantifier) 358.00 (Qualifier)	20	31.20 19.10	15.00 18.00
ZON	12.47	319.1 [M+H] ⁺	283.00 (Quantifier) 187.00 (Qualifier)	20	16.40 25.90	17.00 12.00


LC-MS/MS Information - Typical Chromatograms

Example LC-MS/MS chromatogram for babyfood

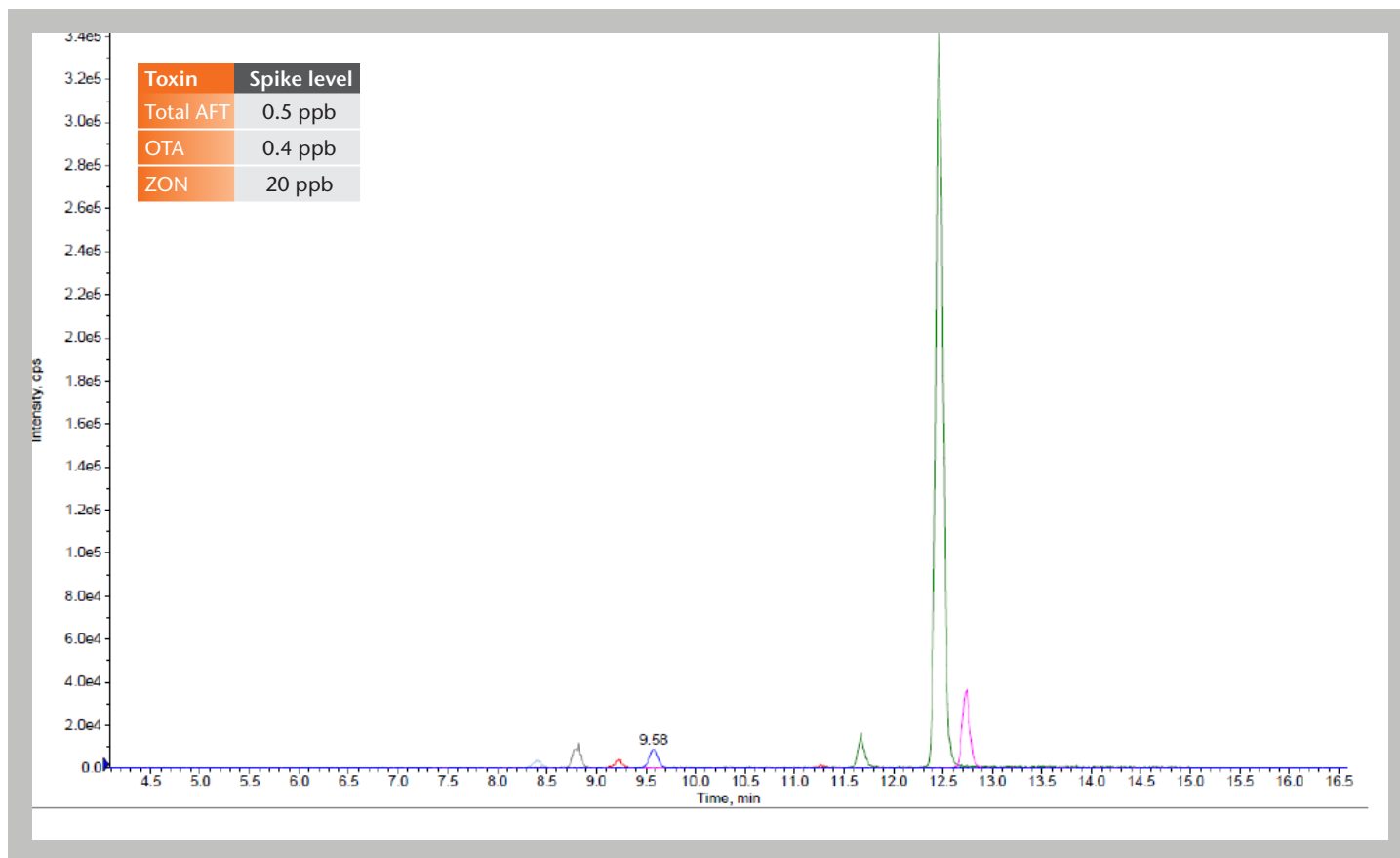

LC-MS/MS Information - Typical Chromatograms

Example LC-MS/MS chromatogram for rye

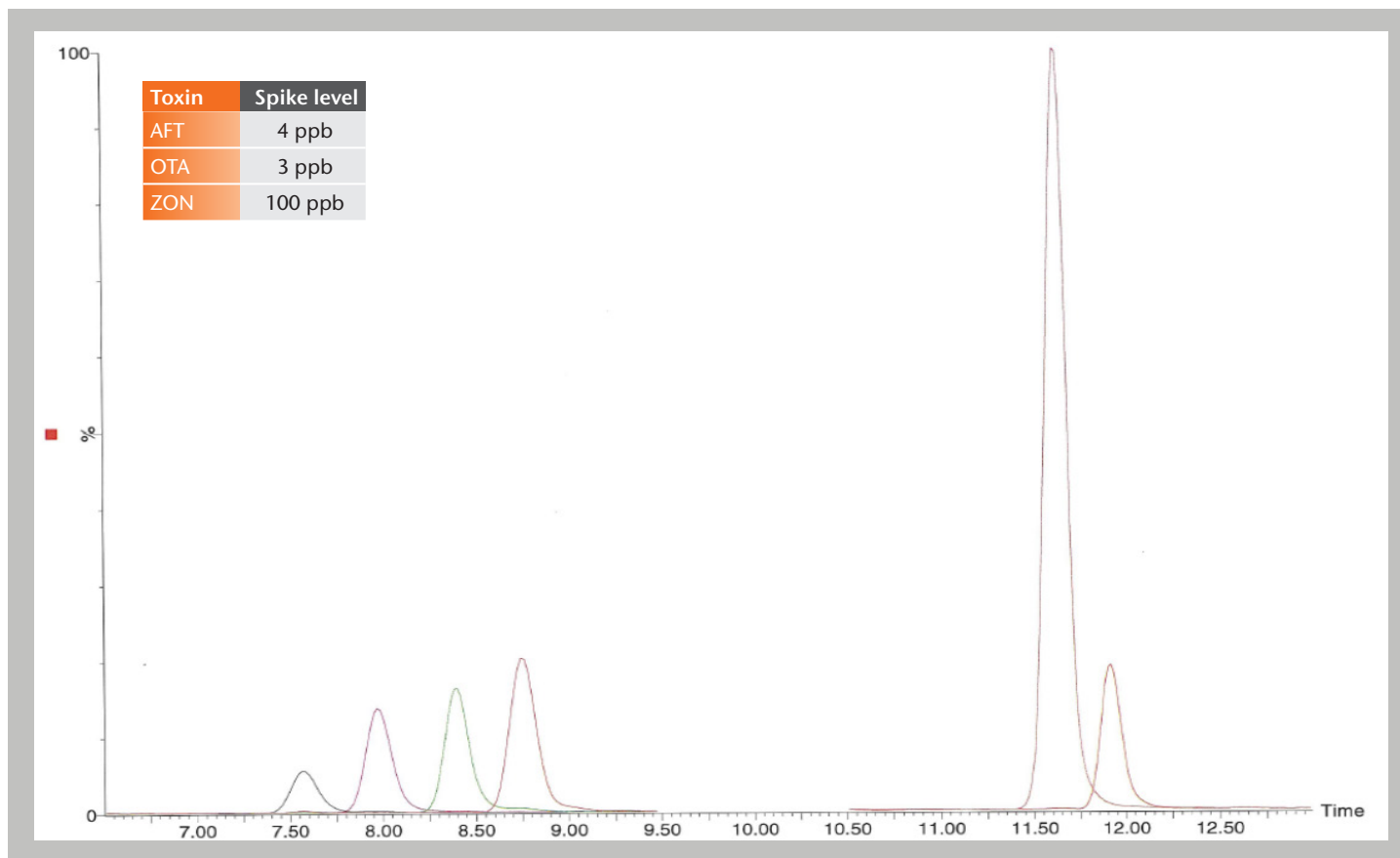

Quality

RBR products are developed, manufactured, tested and dispatched under an ISO 9001 registered Quality Management System, guaranteeing a consistent product, which always meets our performance specifications. Our products have been used in many collaborative studies to develop standard European and International Methods and are widely used by key institutions, food companies and government laboratories. Customer references for RBR products are available on request.

Technical Support

RBR understand that from time to time users of our products may need assistance or advice. Therefore, we are pleased to offer the following services to our customers:

- Analysis of problem samples.
- Application notes for difficult samples.
- References from the RBR library.
- Installation and support of the KOBRA® CELL.
- Advice on detection parameters.
- Advice on preparation and handling of standards.
- Updates on legislation, sampling and other news by e-mail.
- Provision of spiked samples.

Please contact your local R-Biopharm distributor for further information.

Acknowledgement

RBR would like to acknowledge the invaluable help of Peter Mann, Campden BRI, Brewing Division for his assistance during the development of this product.

Warranty

R-Biopharm Rhône Ltd makes no warranty of any kind, express or implied, except that all products made by R-Biopharm Rhône Ltd are made with materials of suitable quality. If any materials are defective, R-Biopharm Rhône Ltd will provide a replacement product. The user assumes all risk and liability resulting from the use of R-Biopharm Rhône Ltd products and procedures. R-Biopharm Rhône Ltd shall not be liable for any damages, including special or consequential damages, loss or expense arising directly or indirectly from the use of R-Biopharm Rhône Ltd products or procedures.

