
CONGEN

SureFast®
Clostridium botulinum
Screening PLUS

Art. No. F5110
100 rxn

User Manual

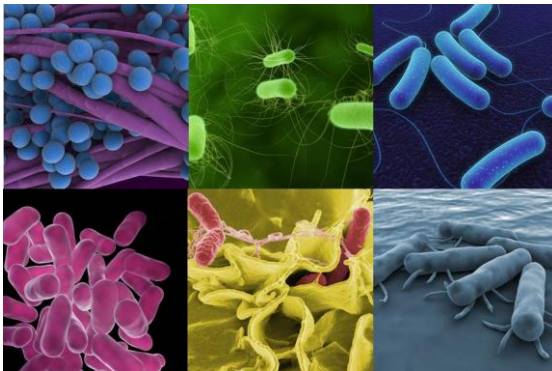April 2024

 **Inhalt**

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Content

1	General Information	8
1.1	Description	8
1.2	Limit of Detection	8
1.3	DNA-preparation	9
1.4	Kit components and storage	9
1.5	Additionally required equipment and materials	9
1.6	Precautions for users	9
1.7	Setup.....	10
1.8	Detection channel Set-up	10
2	Qualitative Analysis	11
2.1	Protocol	11
2.1.1	Preparation of the master-mix	11
2.1.2	Preparation of the real-time PCR-mix	11
2.2	Interpretation of results	12
3	Limitations of the method	12
4	Further Information	12
4.1	Product Information	12
4.2	Technical Support	12

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1 General Information

1.1 Description

The SureFast® Clostridium botulinum Screening PLUS is a real-time PCR for the direct, qualitative detection of a specific Botulinum neurotoxins (BoNT) A, B, E and F DNA sequences of *Clostridium botulinum*, *Clostridium baratii* and *Clostridium butyricum* in enriched food.

Each reaction contains an internal amplification control (IAC). If the DNA contains PCR inhibiting substances, the signal of the amplification control will be affected or the amplification will be suppressed. Examples for PCR inhibiting substances are alcohols (e.g. ethanol, isopropanol), surfactants (e.g. CTAB, SDS, Triton X100) and salts (e.g. sodium chloride). In addition spices, herbs, algae, cocoa and further sample matrices might have PCR inhibiting effects.

The real-time PCR assay can be performed with commonly used real-time PCR instruments, equipped for detection of two fluorescence emissions at the channels FAM and VIC/HEX at the same time. The technical verification of instruments was performed on Roche LightCycler® 480 II, Roche cobas® z 480 Analyzer, Roche LightCycler® 2.0, Qiagen Rotor-Gene Q, Applied Biosystems 7500, Bio-Rad CFX96, R-Biopharm RIDA®CYCLER, LTF MyGo Pro and Agilent Mx3005P.

1.2 Limit of Detection

The SureFast® Clostridium botulinum Screening PLUS real-time PCR has a limit of detection of ≤ 50 DNA copies.

The assay limit of detection depends on sample matrix, processing grade, DNA preparation and DNA content.

The SureFast® PCR systems are very sensitive and therefore even a small amount of target DNA is sufficient for a successful analysis. The concentration of total DNA in the sample does not allow a conclusion on the quantity and quality of the target DNA.

1.3 DNA-preparation

For DNA-preparation the use of SureFast® PREP Bacteria (Art. No. F1021) is recommended.

To assess the process of bacterial growth, it is suggested to compare the samples at the beginning and at the end of the culturing (bacterial growth at Cp difference > 3).

For cultivation the method of the ISO 17919 is recommended.

1.4 Kit components and storage

Kit Code	Reagent	Amount	Lid Color
1	Reaction Mix	2 x 1050 µl	Yellow
2	Taq Polymerase	1 x 80 µl	Dark Red
3	Positive Control	1 x 200 µl	Light Blue

Store all reagents at –28 to - 16°C and protected from light. The Taq Polymerase can be stored at +2 to +8°C for multiple uses on the same day.

Note: The Taq Polymerase may be in a frozen or unfrozen state. This does not affect the quality of the Taq Polymerase or the performance of the real-time PCR.

1.5 Additionally required equipment and materials

- DNA-Extraction kit (e.g. SureFast® PREP Bacteria Art. No. F1021)
- real- time PCR instrument with two detection channels (510 nm and 580 nm)
- real-time PCR consumable (plates, tubes, capillaries, foils, caps)
- pipettes with filter tips
- powder-free disposable gloves
- Vortex mixer
- micro centrifuge with a rotor for the reaction tubes

1.6 Precautions for users

- Extraction, PCR preparation and the PCR run should be separated in different rooms to avoid cross-contaminations.
- This test must only be performed by laboratory personnel trained in molecular biology methods.
- Strictly follow the working instructions.
- When handling samples, wear disposable gloves. After finishing the test, wash your hands.
- Do not smoke, eat or drink in areas where samples or test reagents are being used.
- Food samples and enrichment cultures must be treated as potentially infectious as well as all reagents and materials being exposed to the samples and have to be handled according to the national safety regulations.
- Do not use the kit after the expiration date.

SureFast® Clostridium botulinum Screening PLUS (100 rxn)

Art. No. F5110

April 2024

1.7 Setup

	Blockcycler & R-Biopharm RIDA®CYCLER	Rotorcycler & LightCycler® 480 II & LTF MyGo Pro
Initial Denaturation (HOLD)	1 min, 95°C	1 min, 95°C
Cycles	45	45
Denaturation	15 sec, 95°C	10 sec, 95°C
Annealing/Extension (CYCLE)	30 sec, 60°C	15 sec, 60°C
Temperature Transition Rate/ Ramp Rate	Maximum	Maximum

1.8 Detection channel Set-up

Real-time PCR device	Detection	Detection channel	Quencher	Note
Agilent Mx3005P	<i>Clostridium botulinum</i>	FAM	+	
	IAC	HEX	+	
Applied Biosystems 7500	<i>Clostridium botulinum</i>	FAM	None	Check the passive reference option ROX is none.
	IAC	VIC	None	
Bio-Rad CFX96/Dx/Opus	<i>Clostridium botulinum</i>	FAM	+	
	IAC	VIC/HEX	+	
R-Biopharm RIDA®CYCLER	<i>Clostridium botulinum</i>	green	+	
	IAC	yellow	+	
LTF MyGo Pro	<i>Clostridium botulinum</i>	FAM	+	
	IAC	VIC	+	
Qiagen Rotor-Gene Q	<i>Clostridium botulinum</i>	green	+	
	IAC	yellow	+	
Roche LightCycler® 480 II	<i>Clostridium botulinum</i>	465-510	+	
	IAC	533-580	+	
Roche cobas® z 480 Analyzer	<i>Clostridium botulinum</i>	465-510	+	
	IAC	540-580	+	

2 Qualitative Analysis

2.1 Protocol

2.1.1 Preparation of the master-mix

Calculate the total number of reactions needed (samples and control reactions) for the specific PCR assay. Recommended control reactions for the specific PCR assay: negative control, extraction control, Positive Control. For the analysis of enrichments additional controls are recommended: zero control (sample before enrichment) and medium control.

The reaction mix contains an internal amplification control (inhibition control) per reaction.

Reactions needed for the qualitative *Clostridium botulinum* detection:

3 reactions for controls (1x negative control, 1x extraction control, 1x Positive Control)

For each sample: at least 1 reaction for each sample DNA

Reactions needed for the qualitative *Clostridium botulinum* detection in enrichments:

5 reactions for controls (1x negative control, 1x extraction control, 1x Positive Control, 1x zero control, 1x medium control)

For each sample: at least 1 reaction for each sample DNA

It is also recommended to prepare the master-mix with 10 % additional volume in order to compensate reagent loss. Allow the reagents to thaw, mix and centrifuge before opening and use.

Example for the calculation and preparation of 10 reactions:

Components of the master-mix	Amount per reaction	10 reactions (with 10% excess)
Reaction Mix	19.3 µl	212.3 µl
Taq Polymerase	0.7 µl	7.7 µl
Total volume	20 µl	220 µl

Mix each master-mix well and centrifuge shortly before use.

2.1.2 Preparation of the real-time PCR-mix

- Pipette 20 µl of the master-mix into appropriate tubes/wells.
- Close the negative control (the negative control is ready for PCR without any addition).
- Pipette 5 µl of sample DNA into the designated tubes/wells and close them.
- Pipette 5 µl of Positive Control into the designated tubes/wells and close them.
- Centrifuge all tubes/plates or capillaries shortly at low speed.
- Place tubes/plates into the real-time PCR instrument and start the run according to the setup.

2.2 Interpretation of results

The evaluation has to be made according to the usual analysis program recommended by the real-time PCR instrument manufacturer.

The control reactions have to show the correct results.

Clostridium botulinum DNA is detected in the FAM-channel. In the VIC/HEX-channel the amplification control is detected.

A sample is stated **positive** for *Clostridium botulinum*, if the sample DNA shows amplification in the FAM-channel. High amplicon concentrations can result in a weak or absent signal of the internal amplification control (IAC).

A Cp value for the internal amplification control (IAC) is not needed to obtain a positive result of the Positive Control.

A sample is stated **negative** for *Clostridium botulinum*, if the sample DNA shows no amplification in the FAM-channel and if the internal control (VIC/HEX-channel) of the sample is **positive**.

If the sample DNA in the VIC/HEX-Channel shows **no amplification** or an irregular amplification curve, it contains PCR inhibiting substances. A significant decrease in the fluorescence signal can also show the presence of PCR inhibiting substances. Under these circumstances DNA isolation and purification of the sample need to be improved. Alternatively the DNA can be diluted (recommendation 1:2 in PCR-water) and analysed again for inhibition. Please note that the dilution factor also affects the detection limit of the specific *Clostridium botulinum* PCR assay.

3 Limitations of the method

- The presence of PCR inhibitors may cause invalid results.
- Extremely low levels of target below the limit of detection (LoD) may be detected, but results may not be reproducible.
- A positive test result does not necessarily indicate the presence of a viable organism. It is indicative of the presence of the target DNA (*Clostridium botulinum*).

4 Further Information

4.1 Product Information

- Detailed information about setup of several real-time PCR devices (Download: www.congen.de/en/company/downloads)
- Validation Report upon request

4.2 Technical Support

For further questions please send an e-mail to info@congen.de.