
CONGEN

SureFast[®] Salmonella Species/Enteritidis/Typhimurium 4plex

Art. No. F5166
100 rxn

User Manual

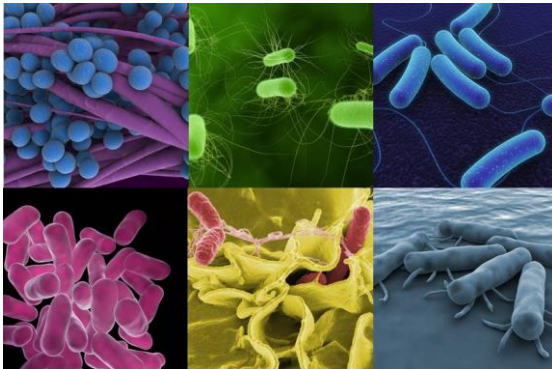January 2025

 **Inhalt**

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Content

1	General Information	10
1.1	Description	10
1.2	Limit of Detection	10
1.3	DNA-preparation	11
1.4	Kit components and storage	11
1.5	Additionally required equipment and materials	11
1.6	Precautions for users	11
1.7	Setup.....	12
1.8	Detection channel Set-up	12
2	Qualitative Analysis	14
2.1	Protocol	14
2.1.1	Preparation of the master-mix	14
2.1.2	Preparation of the real-time PCR-mix	14
2.2	Interpretation of results	15
3	Limitations of the method	16
4	Further Information	16
4.1	Product Information	16
4.2	Technical Support	16

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1 General Information

1.1 Description

The SureFast® Salmonella Species/Enteritidis/Typhimurium 4plex is a real-time PCR for the direct, qualitative detection and differentiation of specific DNA sequences of *Salmonella* spp., *Salmonella* Enteritidis and *Salmonella* Typhimurium.

Each reaction contains an internal amplification control (IAC). If the DNA contains PCR inhibiting substances, the signal of the amplification control will be affected or the amplification will be suppressed. Examples for PCR inhibiting substances are alcohols (e.g. ethanol, isopropanol), surfactants (e.g. CTAB, SDS, Triton X100) and salts (e.g. sodium chloride). In addition spices, herbs, algae, cocoa and further sample matrices might have PCR inhibiting effects.

The real-time PCR assay can be performed with commonly used real-time PCR instruments, equipped for detection of four fluorescence emissions at the channels FAM, VIC/HEX, ROX and Cy5 at the same time. The technical verification of instruments was performed on Agilent AriaDx, Agilent Mx3005P, Applied Biosystems 7500, Bio-Rad CFX96, Bio Molecular Systems MIC, Qiagen Rotor-Gene Q, Roche cobas® z 480 Analyzer and Roche LightCycler® 480 II.

1.2 Limit of Detection

The SureFast® Salmonella Species/Enteritidis/Typhimurium 4plex real-time PCR has a limit of detection of ≤ 5 DNA copies.

The assay limit of detection depends on sample matrix, processing grade, DNA preparation and DNA content.

The SureFast® PCR systems are very sensitive and therefore even a small amount of target DNA is sufficient for a successful analysis. The concentration of total DNA in the sample does not allow a conclusion on the quantity and quality of the target DNA.

SureFast® Salmonella Species/Enteritidis/Typhimurium 4plex (100 rxn)

Art. No. F5166

January 2025

1.3 DNA-preparation

For DNA-preparation the use of SureFast® PREP Bacteria (Art. No. F1021) or SureFast® Speed PREP (Art. No. F1054) is recommended.

To assess the process of bacterial growth, it is suggested to compare the samples at the beginning and at the end of the culturing (bacterial growth at Cp difference > 3).

For cultivation the method of the ISO 10135 is recommended.

1.4 Kit components and storage

Kit Code	Reagent	Amount	Lid Color
1	Reaction Mix	2 x 1050 µl	Yellow
2	Taq Polymerase	1 x 80 µl	Dark Red
3	Positive Control	1 x 200 µl	Light Blue

Store all reagents at -28 to -16 °C and protected from light. The Taq Polymerase can be stored at +2 to +8 °C for multiple uses on the same day.

Note: The Taq Polymerase may be in a frozen or unfrozen state. This does not affect the quality of the Taq Polymerase or the performance of the real-time PCR.

1.5 Additionally required equipment and materials

- DNA-Extraction kit
(e.g. SureFast® PREP Bacteria Art. No. F1021 / SureFast® Speed PREP Art. No. F1054)
- real-time PCR instrument with four detection channels (510 nm, 580 nm, 610 nm and 660 nm)
- real-time PCR consumable (plates, tubes, foils, caps)
- pipettes with filter tips
- powder-free disposable gloves
- Vortex mixer
- microcentrifuge with a rotor for the reaction tubes

1.6 Precautions for users

- Extraction, PCR preparation and the PCR run should be separated in different rooms to avoid cross-contaminations.
- This test must only be performed by laboratory personnel trained in molecular biology methods.
- Strictly follow the working instructions.
- When handling samples, wear disposable gloves. After finishing the test, wash your hands.
- Do not smoke, eat or drink in areas where samples or test reagents are being used.
- Food samples and enrichment cultures must be treated as potentially infectious as well as all reagents and materials being exposed to the samples and have to be handled according to the national safety regulations.
- Do not use the kit after the expiration date.

SureFast® Salmonella Species/Enteritidis/Typhimurium 4plex (100 rxn)

Art. No. F5166

January 2025

1.7 Setup

	Blockcyler & R-Biopharm RIDA®CYCLER	Rotorcyler & LightCycler® 480 II
Initial Denaturation (HOLD)	1 min, 95°C	1 min, 95°C
Cycles	45	45
Denaturation	15 sec, 95°C	10 sec, 95°C
Annealing/Extension (CYCLE)	30 sec, 60°C	15 sec, 60°C
Temperature Transition Rate/ Ramp Rate	Maximum	Maximum

1.8 Detection channel Set-up

Real-time PCR device	Detection	Detection channel	Quencher	Note
Agilent Mx3005P	<i>Salmonella</i> spp.	FAM	+	
	IAC	HEX	+	
	<i>Salmonella</i> Enteritidis	ROX	+	
	<i>Salmonella</i> Typhimurium	Cy5	+	
Agilent AriaDx / Mx	<i>Salmonella</i> spp.	FAM	+	
	IAC	HEX	+	
	<i>Salmonella</i> Enteritidis	ROX	+	
	<i>Salmonella</i> Typhimurium	Cy5	+	
Applied Biosystems 7500	<i>Salmonella</i> spp.	FAM	None	Check the passive reference option ROX is none.
	IAC	VIC	None	
	<i>Salmonella</i> Enteritidis	ROX	None	
	<i>Salmonella</i> Typhimurium	Cy5	None	
Bio-Rad CFX96	<i>Salmonella</i> spp.	FAM	+	
	IAC	VIC/HEX	+	
	<i>Salmonella</i> Enteritidis	ROX	+	
	<i>Salmonella</i> Typhimurium	Cy5	+	
R-Biopharm RIDA®CYCLER/ Bio Molecular Systems MIC	<i>Salmonella</i> spp.	green	+	Ignore cycles before , if there is a significant deviation in the baseline at the start of the run. Please see page 45 of the cyler operating instructions, section 12.1.2 Cycling analysis parameter.
	IAC	yellow	+	
	<i>Salmonella</i> Enteritidis	orange	+	
	<i>Salmonella</i> Typhimurium	red	+	

**SureFast® Salmonella
Species/Enteritidis/Typhimurium 4plex
(100 rxn)**

Art. No. F5166

January 2025

Real-time PCR device	Detection	Detection channel	Quencher	Note
Qiagen Rotor-Gene Q	<i>Salmonella</i> spp.	green	+	Note: Please use only 0.1 ml reaction tubes. The gain settings must be set to 5 (factory default) for all channels.
	IAC	yellow	+	
	<i>Salmonella</i> Enteritidis	orange	+	
	<i>Salmonella</i> Typhimurium	red	+	
Roche LightCycler® 480 II	<i>Salmonella</i> spp.	465-510	+	The SureCC Color Compensation Kit I (Art. No. F4009) is required.
	IAC	533-580	+	
	<i>Salmonella</i> Enteritidis	533-610	+	
	<i>Salmonella</i> Typhimurium	618-660	+	
Roche cobas® z 480 Analyzer	<i>Salmonella</i> spp.	465-510	+	The SureCC Color Compensation Kit I (Art. No. F4009) is required.
	IAC	540-580	+	
	<i>Salmonella</i> Enteritidis	540-610	+	
	<i>Salmonella</i> Typhimurium	610-670	+	

2 Qualitative Analysis

2.1 Protocol

2.1.1 Preparation of the master-mix

Calculate the total number of reactions needed (samples and control reactions) for the specific PCR assay as well as for the inhibition control.

Recommended control reactions for the specific PCR assay: negative control, extraction control and Positive control.

For the analysis of enrichments additional controls are recommended: zero control (sample before enrichment) and medium control.

The reaction mix contains an internal amplification control (Inhibition control) per reaction.

Reactions needed for the qualitative *Salmonella* detection:

3 reactions for controls (1x negative control, 1x extraction control, 1x Positive Control)

For each sample: at least 1 reaction for each sample DNA

Reactions needed for the qualitative *Salmonella* detection in enrichments:

5 reactions for controls (1x negative control, 1x extraction control, 1x Positive Control, 1x zero control, 1x medium control)

For each sample: at least 1 reaction for each sample DNA

It is also recommended to prepare the master-mix with 10% additional volume in order to compensate reagent loss. Allow the reagents to thaw, mix and centrifuge before opening and use.

Example for the calculation and preparation of 10 reactions:

Components of the master-mix	Amount per reaction	10 reactions (with 10% excess)
Reaction Mix	19.3 µl	212.3 µl
Taq Polymerase	0.7 µl	7.7 µl
Total volume	20 µl	220 µl

Mix each master-mix well and centrifuge shortly before use.

2.1.2 Preparation of the real-time PCR-mix

- Pipette 20 µl of the master-mix into appropriate tubes/wells.
- Close the negative control (the negative control is ready for PCR without any addition).
- Pipette 5 µl of sample DNA into the designated tubes/wells and close them.
- Pipette 5 µl of Positive Control into the designated tubes/wells and close them.
- Centrifuge all tubes/plates shortly at low speed.
- Place tubes/plates into the real-time PCR instrument and start the run according to the setup.

2.2 Interpretation of results

The evaluation has to be made according to the usual analysis program recommended by the real-time PCR instrument manufacturer.

The control reactions have to show the correct results.

Salmonella spp. DNA is detected in the FAM-channel, *Salmonella* Enteritidis DNA is detected in the ROX-channel and *Salmonella* Typhimurium DNA is detected in the Cy5-channel (see table). In the VIC/HEX-channel the internal amplification control (IAC) is detected.

A sample is stated **positive** for the respective parameter, if the sample DNA shows amplification in the respective channel. High amplicon concentrations can result in a weak or absent signal of the IAC.

A Cp value for the IAC is not needed to obtain a positive result of the Positive Control.

A sample is stated **negative** for the respective parameter, if the sample DNA shows no amplification in the respective channel and if the internal control (VIC/HEX-channel) of the sample is **positive** with a shift in Cp-Value ≤ 2 compared to the negative control.

If the sample DNA in the VIC/HEX-Channel shows **no amplification** or a shift in Cp-value > 2 compared to the negative control, it contains PCR inhibiting substances. A significant decrease in the fluorescence signal can also show the presence of PCR inhibiting substances. Under these circumstances DNA isolation and purification of the sample need to be improved. Alternatively the DNA can be diluted (recommendation 1:2 in PCR-water) and analysed again for inhibition. Please note that the dilution factor also affects the detection limit of the specific *Salmonella* spp., *Salmonella* Enteritidis or *Salmonella* Typhimurium PCR assay.

Result in the respective channel				Interpretation
FAM channel <i>Salmonella</i> spp.	ROX channel <i>Salmonella</i> Enteritidis	Cy5 channel <i>Salmonella</i> Typhimurium	VIC/HEX channel IAC	
positive	negative	negative	positive	<i>Salmonella</i> spp. DNA detected
positive	positive	negative	positive	<i>Salmonella</i> Enteritidis DNA detected
positive	negative	positive	positive	<i>Salmonella</i> Typhimurium DNA detected
negative	negative	negative	negative	invalid

Note: The results displayed in the table above represent merely an example. Additional combinations are also possible.

3 Limitations of the method

- The presence of PCR inhibitors may cause invalid results.
- Extremely low levels of target below the limit of detection (LoD) may be detected, but results may not be reproducible.
- A positive test result does not necessarily indicate the presence of viable organism. It is indicative for the presence of the target DNA.
- Cross reactivity was observed with DNA extracts from *Salmonella* Amsterdam, *Salmonella* Blegdam, *Salmonella* Kiel, *Salmonella* Moscow and *Salmonella* Nitra in the detection system for *Salmonella* Enteritidis.

4 Further Information

4.1 Product Information

- Detailed information about setup of several real-time PCR devices
(Download: www.congen.de/en/downloads)
- Product-related documents (Download: www.congen.de/en/eifu/)
- Validation Report upon request

4.2 Technical Support

For further questions please send an e-mail to info@congen.de.