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QuickGEN PCR Kit

P1 Screening and Hop resistance

1. Examples of successful detected microorganisms

- *Lactobacillus* spp.
 - L. acetotolerans*
 - L. acidophilus*
 - L. amylolyticus*
 - L. amylovorus*
 - L. backii*
 - L. brevis*
 - L. brevisimilis*
 - L. buchneri*
 - L. casei* / *paracasei*
 - L. collinoides* / *paracasei*
 - L. coryniformis* ssp. *coryniformis*
 - L. coryniformis* ssp. *torquens*
 - L. curvatus*
 - L. delbrückii* ssp. *delbrückii*/ *lactis*
 - L. harbinensis*
 - L. helveticus*
 - L. johnsonii*
 - L. lindneri*
 - L. malefermentans*
 - L. parabrevis*
 - L. parabuchneri* (*frigidus*)
 - L. paucivorans*
 - L. pentosus*
 - L. perolens*
 - L. plantarum*
 - L. paraplantarum*
 - L. rossiae*
- *Lactococcus lactis*
- *Megasphaera* spp.
- *Micrococcus kristinae*
- *Pectinatus* spp.
 - Pectinatus cerevisiophilus*
 - Pectinatus frisingensis*
 - Pectinatus haikarae*
 - Pectinatus* ssp.
- *Pediococcus* spp.
 - Pediococcus acidilactici*
 - Pediococcus damnosus*
 - Pediococcus dextrinicus*
 - Pediococcus inopinatus*
 - Pediococcus pentosaceus*
 - Pediococcus parvulus*
 - Pediococcus claussenii*

2. Intended use

Qualitative detection of bacterial contaminations and the hop resistance genes *horA* and *horC* as genetic markers for beer-spoilage microorganisms in beer and beer mixes.

3. Test principle

The TaqMan® real-time PCR is based on hot-start-PCR and sequence-specific dual labelled probes (FAM/DQ; HEX/DQ; ROX/DQ) which, when accurately hybridised, emit a measurable fluorescent signal of a defined wavelength in the extension phase. The increase of signal is continuously measured in a real-time PCR detection instrument. To avoid false negative PCR-results an Inhibition Control (CY5/DQ) is amplified together in one reaction vessel with the specific sequence.

4. Kit contents

The kit contains sufficient reagents for 48 reactions:

1 x Premix	white cap
48 x Dye Strips (freeze-dried, incl. IC-DNA)	tube strips
6 x Cap Strips	
1 x ddH ₂ O	colourless cap
1 x Control-DNA (freeze-dried)	yellow cap

5. Storage conditions

The Control-DNA is freeze-dried and has to be solved in ddH₂O prior to use (see 7.1).

Do **not** freeze the lyophilized Dye Strips and lyophilized Control-DNA.

The PCR-reagents should be stored at 2 – 8 °C (35 – 46 °F).

Keep Premix for storage at - 20 °C (- 4 °F) **after arrival**. Avoid loss of sensitivity by repeating freezing and thawing more than 3 times. For irregular use aliquot the Premix.

The Dye Strips contain the fluorescent labelled probes and should be handled light protected. All reagents are stable for 12 months, if they are stored correctly.

6. Materials required, but not provided

6.1. Instruments

Real-time PCR machine for high profile PCR tubes

Standard benchtop mini-centrifuge (1.5 – 2.0 mL vials)

Centrifuge for strips

Pipettes

“Vortex”

6.2. Reagents and plastic ware

sterile ddH₂O

sterile filter tips

7. PCR

7.1. PCR-Setup

When using the kit for the first time, the freeze-dried Control-DNA has to be shortly centrifuged and carefully resolved:

- add 55 µL sterile ddH₂O to the freeze-dried Control-DNA
- after 15 minutes mix well

Before every use thoroughly mix all PCR-components and centrifuge briefly.

Remove the film from the needed tube strips and pipette the PCR-components. After pipetting close the tube strips with the provided cap strips.

PCR-reaction Setup:

PCR-components	volume (µL)
Premix	15.0
Sample-DNA	5.0*
Total volume	20.0

* if using the Simplex Easy DNA- or QuickGEN Yeast Sample Preparation Kit:
add 2.5 µL DNA and 2.5 µL sterile ddH₂O

1. Pipette 15 µL of the Premix into each PCR-tube, making sure that, prior to the first filling, the tip has been moistened.
2. Add 5 µL sample DNA, add 5 µL of the Control-DNA for the positive control reaction, add 5 µL of the extraction control and 5 µL of ddH₂O for the negative control* reaction. Use a fresh tip with each DNA filling.
3. Close the tubes immediately and centrifuge them shortly.
4. Place the tubes in the PCR-machine and start run.

Very important: * Please fill up the negative control with 5 µL ddH₂O to avoid unspecific amplification.

Work swiftly to avoid warming up and keep away from light.

7.2 PCR-Program

The Premix contains no ROX. This must be considered according to the settings of the real-time machine. **ABI 7500:** Under “Assign Targets and Samples” in ‘Select the dye to use as the passive reference’ choose „none“.

Mx3500P: Under Instruments in ‘Filter Set Gain Settings’ change the Filter set gain multiplier to CY5 x1, ROX x1, HEX-JOE x4 und FAM x8.

Step	Time	Temp.	
Initial denaturation	15 min	95 °C	
Cycling Denaturation	10 sec	95 °C	Cycle 35 x
Cycling Annealing/ Elongation	20 sec	60 °C	

8. Evaluation

The evaluation has to be made according to the data analysis program recommended by the real-time instrument manufacturer.

Lactobacillus/ Pediococcus-DNA: FAM-channel

Megasphaera/ Pectinatus-DNA: ROX-channel

horA and horC: HEX-channel

Inhibition Control-DNA: CY5-channel


A sample is ***Lactobacillus/ Pediococcus* positive**, if there is a detectable fluorescence increase in the **FAM-channel** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls it has to be positive.

A sample is ***Megasphaera/ Pectinatus* positive**, if there is a detectable fluorescence increase in the **ROX-channel** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls it has to be positive.

A sample is **horA and/or horC positive**, if there is a detectable fluorescence increase in the **HEX-channel** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls has to be positive.

A sample is **negative**, if there is no detectable fluorescence increase in the different channels and the positive controls have a positive fluorescence signal. The negative controls show no amplification. The Inhibition Control in the CY5-channel has to be positive in the sample and in the negative controls, a false negative result due to inhibitory effects is then excluded.

analysis flowchart

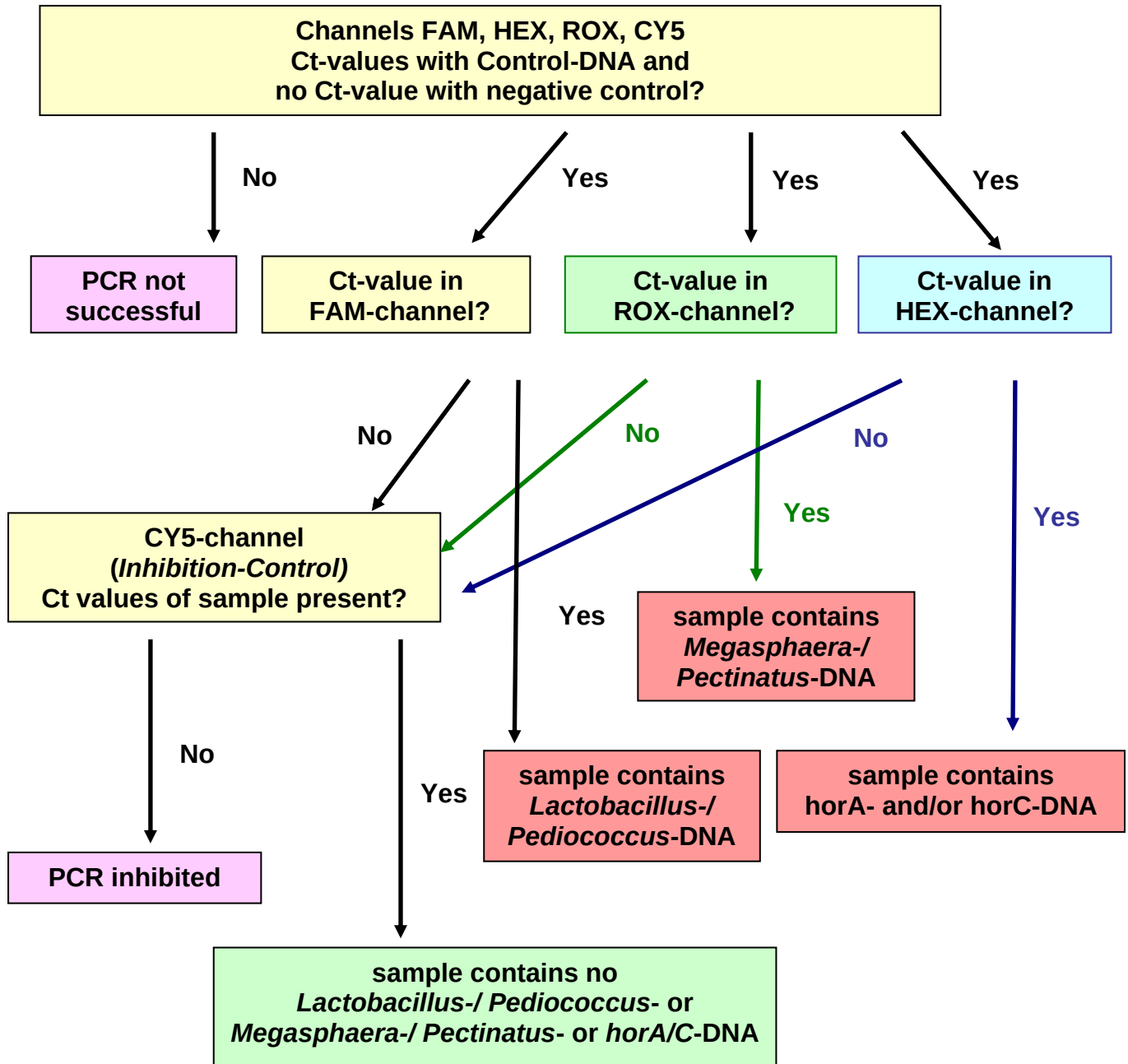Note:

The polymerase-chain reaction (PCR) is protected by patents and requires a licence from Hoffmann-LaRoche Inc..

The provided product does not authorise the purchaser for the commercial use of this method.

GEN-IAL makes no warranty of any kind, either expressed or implied, except that the materials from which its products are made are of standard quality. If any materials are defective, GEN-IAL will provide a replacement product. There is no warranty of merchantability of this product, or of the fitness of the product for any purpose. GEN-IAL shall not be liable for any damages, including special or consequential damage, or expense arising directly or indirectly from the use of this product.