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QuickGEN PCR Kit

P1 Screening and Hop resistance

1. Examples of successful detected microorganisms

- *Lactobacillus* spp.
 - L. acetotolerans*
 - L. acidophilus*
 - L. amylolyticus*
 - L. amylovorus*
 - L. backii*
 - L. brevis*
 - L. brevisimilis*
 - L. buchneri*
 - L. casei* / *paracasei*
 - L. collinoides* / *paracasei*
 - L. coryniformis* ssp. *coryniformis*
 - L. coryniformis* ssp. *torquens*
 - L. curvatus*
 - L. delbrückii* ssp. *delbrückii*/ *lactis*
 - L. harbinensis*
 - L. helveticus*
 - L. johnsonii*
 - L. lindneri*
 - L. malefermentans*
 - L. parabrevis*
 - L. parabuchneri* (*frigidus*)
 - L. paucivorans*
 - L. pentosus*
 - L. perolens*
 - L. plantarum*
 - L. paraplantarum*
 - L. rossiae*
- *Lactococcus lactis*
- *Megasphaera* spp.
- *Micrococcus kristinae*
- *Pectinatus* spp.
 - Pectinatus cerevisiiphilus*
 - Pectinatus frisingensis*
 - Pectinatus haikarae*
 - Pectinatus* ssp.
- *Pediococcus* spp.
 - Pediococcus acidilactici*
 - Pediococcus damnosus*
 - Pediococcus dextrinicus*
 - Pediococcus inopinatus*
 - Pediococcus pentosaceus*
 - Pediococcus parvulus*
 - Pediococcus claussenii*

2. Intended use

Qualitative detection of bacterial contaminations and the hop resistance genes *horA* and *horC* as genetic markers for beer-spoilage microorganisms in beer and beer mixes.

3. Test principle

The TaqMan® real-time PCR is based on hot-start-PCR and sequence-specific dual labelled probes (FAM/DQ; HEX/DQ; ROX/DQ) which, when accurately hybridised, emit a measurable fluorescent signal of a defined wavelength in the extension phase. The increase of signal is continuously measured in a real-time PCR detection instrument. To avoid false negative PCR-results an Inhibition Control (CY5/DQ) is amplified together in one reaction vessel with the specific sequence.

4. Kit contents

The kit contains sufficient reagents for 50 reactions:

1 x Premix	white cap
1 x Dye Mix (freeze-dried, incl. IC-DNA)	dark vial, red cap
1 x ddH ₂ O	colourless cap
1 x Control-DNA (freeze-dried)	yellow cap

5. Storage conditions

The Dye Mix and the Control-DNA are freeze-dried, they have to be solved in ddH₂O prior to use (see 7.1).

Do **not** freeze the lyophilized Dye Mix and lyophilized Control-DNA.

The PCR reagents should be stored at 2 – 8 °C (35 – 46 °F).

Keep Premix for storage at - 20 °C (- 4 °F) **after arrival**. Avoid loss of sensitivity by repeating freezing and thawing more than 3 times. For irregular use aliquot the Premix.

The Dye Mix contains the fluorescent labelled probes and should be handled light protected.

All reagents are stable for 12 months, if they are stored correctly

6. Materials required, but not provided

6.1. Instruments

Real-time PCR instrument with channels FAM, HEX, ROX and CY5

Centrifuge for 1.5 – 2.0 mL reaction tubes

Centrifuge for multiwell plates or strips

Pipettes

“Vortex”

6.2. Reagents and plastic ware

sterile ddH₂O

sterile reaction vessels 1.5 – 2.0 mL

sterile plates or optical tubes (plastic)

sterile filter tips

optional: Colour Compensation Set for LC480 (GEN-IAL[®], REF: Q800)

LC480: Use a Colour compensation prior to first use of the kit

7. PCR

7.1. PCR-Setup

When using the kit for the first time, the freeze-dried kit components have to be shortly centrifuged and carefully resolved:

- add 80 µL sterile ddH₂O to freeze-dried Dye Mix
- add 55 µL sterile ddH₂O to the freeze-dried Control-DNA
- after 15 minutes mix well

Before every use thoroughly mix all PCR-components and centrifuge briefly.

PCR-reaction Setup:

PCR-Components	volume (µL)
Premix	13.5
Dye Mix	1.5
Sample-DNA	5.0*
Total volume	20.0

* if using the Simplex Easy DNA- or QuickGEN Sample Preparation Centrifugation Kits:
add 2.5µL DNA and 2.5 µL sterile ddH₂O

1. Prepare a mastermix by mixing Premix and Dye Mix
2. Multiply said volumes with the number of PCR preparations including controls (positive control, negative control, extraction control), taking into account pipette reserves of approximately 5-10 %.
3. Divide 15 µL of the PCR-mastermix among the individual reaction vessels, making sure that, prior to the first filling, the tip of the pipette has been moistened.
4. Add 5 µL sample DNA, add 5 µL of the Control-DNA for the PCR positive control, add 5 µL of the extraction control and 5 µL of ddH₂O for the negative control* reaction. Use a fresh tip with each DNA filling.
5. Close the tubes immediately and centrifuge them shortly.
6. Place the tubes in the PCR-machine and start run.

Very important: * Please fill up the negative control with 5 µL ddH₂O to avoid unspecific amplification.

Work swiftly to avoid warming up and keep away from light.

7.2 PCR-Program

7.2.1 PCR-Program LC480

1. Click the button **tool** on the right side in the window **LightCycler 480 Software release 1.5.0. SP1**
2. Click the button **Detection formats** at the left side of the menu bar
3. Click **New** in the window **Detection formats** and name the experiment
4. Open the window **Filter Combination Selection** and choose the following filter combinations: 465-510 / 533-580 / 533-610 / 618-660.
5. Open the window **Selected Filter Combination** list and add the following amounts

Excitation Filter	Emission Filter	Name	Melt factor	Quant factor	Max. Integ. Time
465	510	465-510 FAM	1	10	2
533	580	533-580 HEX	1	10	2
533	610	533-610 610	1	10	2
618	660	618-660 660	1	10	2

6. **Close** the window
7. Click the button **New Experiment** on the right side of the menu bar
8. From the pull-down menu **Detection formats** choose the defined experiment, click the button **Customize** and check the detection formats. All of them have to be activated
9. Click the button **ok**
10. Define the following program:

1. Program Name: **Heat**

Cycles **1** Analysis Mode **None**

Target (°C)	Acquisition Mode	Hold (hh:mm:ss)	Ramp Rate (°C/s)	Acquisitions (per °C)	Sec Target (°C)	Step size (°C)	Step Delay (cycles)
95	None	00:15:00	4.40		0	0	0

2. Program Name: **Ampli**

Cycles **30*-35** Analysis Mode **Quantification**

Target (°C)	Acquisition Mode	Hold (hh:mm:ss)	Ramp Rate (°C/s)	Acquisitions (per °C)	Sec Target (°C)	Step size (°C)	Step Delay (cycles)
95	None	00:00:10	4.40		0	0	0
60	Single	00:00:20	2.20		0	0	0

*** for preenriched samples program 30 cycles**

3. Program name : **Cool**

Cycles **1** Analysis Mode **None**

Target (°C)	Acquisition Mode	Hold (hh:mm:ss)	Ramp Rate (°C/s)	Acquisitions (per °C)	Sec Target (°C)	Step size (°C)	Step Delay (cycles)
40	None	00:00:20	2.20		0	0	0

Optional: Saving the program as **run template**. The template button allows to select and apply a template to the currently open object and to save the currently open object as a template. Click the clamp beside the button **Apply Template**. Click **Save as template** and save the file in **templates**.

11. Click the button **Subset editor** on the left side

12. Click the button **+** and **New Subset 1** appears

13. Mark the wells in the **New Subset 1 Settings** window

14. Click the button **Apply**

15. Click the button **Sample editor** on the left side

16. **Very important:** Activate In the window **Step 1 Select Workflow: Abs.Quant**

17. Choose the subset **New Subset** in the window **Step 2 Select Samples**

18. Define your probes in the **Sample table**

19. Switch to the button **Experiment** and start the run with the button **start run**

7.2.2 PCR program for other real-time machines

The Premix contains no ROX. This must be considered according to the settings of the real-time machine. **ABI 7500:** Under “Assign Targets and Samples” in ‘Select the dye to use as the passive reference’ choose „none“.

Step	Time	Temp.	
Initial denaturation	15 min	95 °C	
Cycling Denaturation	10 sec	95 °C	Cycle 30 ¹ - 35 x
Cycling Annealing/ Elongation	20-30 ² sec	60 °C	

¹ for preenriched samples program 30 cycles

² please program 20 sec, if the real-time-PCR device allows it

8. Evaluation

The evaluation has to be made according to the data analysis program recommended by the real-time instrument manufacturer.

For LC480: For analysis activate Colour Compensation

Lactobacillus/ Pediococcus-DNA: FAM-channel (LC480: 465-510 nm)

Megasphaera/ Pectinatus-DNA: ROX-channel (LC480: 533-610 nm)

horA and horC: HEX-channel (LC480: 533-580 nm)

Inhibition Control-DNA: CY5-channel (LC480: 618-660 nm)

A sample is ***Lactobacillus/ Pediococcus* positive**, if there is a detectable fluorescence increase in the **FAM-channel (LC480: 465-510 nm)** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel (LC480: 618-660 nm) may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls it has to be positive.


A sample is ***Megasphaera/ Pectinatus* positive**, if there is a detectable fluorescence increase in the **ROX-channel (LC480: 533-610 nm)** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel (LC480: 618-660 nm) may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls it has to be positive.

A sample is **horA and/or horC positive**, if there is a detectable fluorescence increase in the **HEX-channel (LC480: 533-580 nm)** and the negative controls show no amplification. The positive controls should have a positive fluorescence signal. The Inhibition Control in the CY5-channel (LC480: 618-660 nm) may be positive or negative (depending on the amount of DNA or inhibitors in the sample reaction). In the negative controls has to be positive.

A sample is **negative**, if there is no detectable fluorescence increase in the different channels and the positive controls have a positive fluorescence signal. The negative controls show no amplification. The Inhibition Control in the CY5-channel (LC480: 618-660 nm) has to be positive in the sample and in the negative controls, a false negative result due to inhibitory effects is then excluded.

analysis flowchart

LC480: analysis after activating Colour Compensation

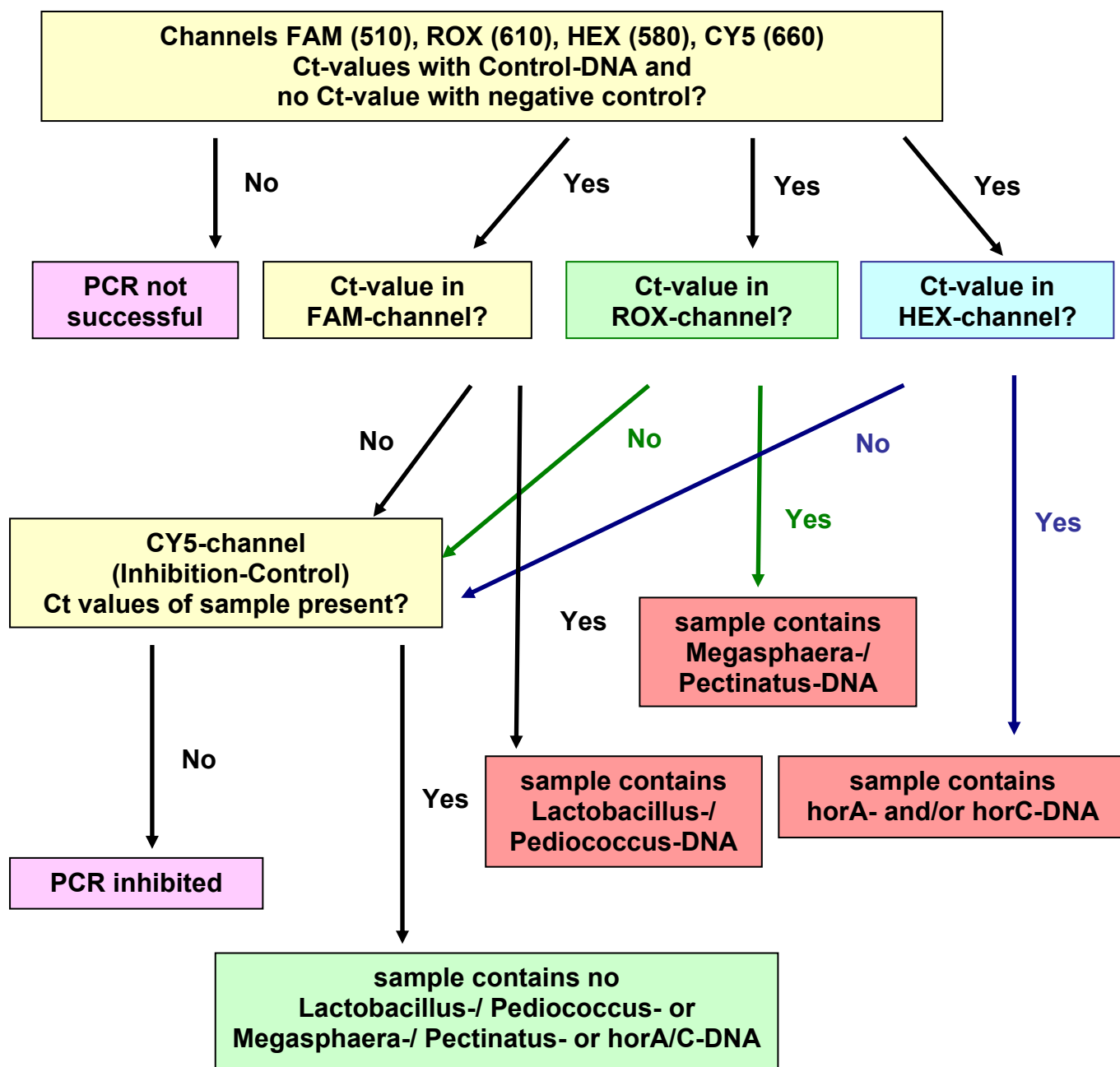Note:

The polymerase-chain reaction (PCR) is protected by patents and requires a licence from Hoffmann-LaRoche Inc..

The provided product does not authorise the purchaser for the commercial use of this method.

GEN-IAL makes no warranty of any kind, either expressed or implied, except that the materials from which its products are made are of standard quality. If any materials are defective, GEN-IAL will provide a replacement product. There is no warranty of merchantability of this product, or of the fitness of the product for any purpose. GEN-IAL shall not be liable for any damages, including special or consequential damage, or expense arising directly or indirectly from the use of this product.