# RIDA® QUICK Aflatoxin RQS ECO

**REF** R5209

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# RIDA<sup>®</sup>QUICK Aflatoxin RQS ECO

## Brief information

RIDA<sup>®</sup>QUICK Aflatoxin RQS ECO (Art. No. R5209) is an immunochromatographic test for the quantitative determination of aflatoxin in corn. The test uses an aqueous extraction method.

The test kit contains 20 test strips, each of which can be used for one analysis. All reagents required for the assay are contained in the test kit. Results are evaluated with

- the RIDA<sup>®</sup>SMART APP software (Art. No. ZRSAM) and an approved smartphone or
- the RIDA<sup>®</sup>SMART APP software (Art. No. ZRSAM) installed on an Android smartphone and the RIDA<sup>®</sup>SMART BOX (Art. No. ZRSA-SB).

Sample preparation:	homogenization and extraction
Time requirement:	sample preparation (for 10 samples) ..... approx. 10 min test implementation (incubation time) ..... 5 min
Limit of detection:	..... < 2 µg/kg (ppb)
Limit of quantification:	..... 2.89 µg/kg
Measuring range:	2 - 75 µg/kg and 50 - 300 µg/kg Quantitative evaluation with the RIDA <sup>®</sup> SMART APP software
Specificity:	The RIDA <sup>®</sup> QUICK Aflatoxin RQS ECO test detects aflatoxin in corn.

Further information is contained in the validation report.

The specificity of the RIDA<sup>®</sup>QUICK Aflatoxin RQS ECO test was established by determining the cross-reactivity to the corresponding substances in the buffer system. In samples, the specificity may deviate from those determined in the buffer system due to matrix effects. Prior to the analysis of cross-reactive substances, the user has to determine the Limit of Detection and the Recovery for the substance in the respective sample matrix. The test cannot discriminate between analytes and cross-reactive substances.

## Related products and accessories for Aflatoxin determination

RIDASCREEN® Aflatoxin Total	(Art. No. R4701)
RIDASCREEN® Aflatoxin B1 30/15	(Art. No. R1211)
RIDASCREEN® Aflatoxin M1	(Art. No. R1121)
RIDASCREEN®FAST Aflatoxin	(Art. No. R5202)
RIDASCREEN®FAST Aflatoxin SC	(Art. No. R9002)
RIDA®QUICK Aflatoxin RQS	(Art. No. R5208)
RIDA® Aflatoxin column	(Art. No. R5001 / R5002)
RIDA®SMART APP Software	(Art. No. ZRSAM)
RIDA®SMART BOX	(Art. No. ZRSA-SB)
Trilogy® certified reference material and standards for Aflatoxin (ISO 17034)	

### 1. Intended use

RIDA®QUICK Aflatoxin RQS ECO is an immunochromatographic test in strip format for the quantitative determination of aflatoxin in corn. The test uses an aqueous extraction method of aflatoxin.

### 2. General information

Aflatoxins are secondary metabolites of the fungi species *Aspergillus flavus*, *parasiticus* and *nomius*. These fungal species are found in humid, tropical areas and can contaminate plant-based food and feed. This contamination usually takes place in the countries of cultivation. Aflatoxins belong to the strongest natural occurring carcinogenic substances.

### 3. Test principle

The immunochromatographic test in the form of a test strip is based on an antigen-antibody reaction. A specific anti-aflatoxin antibody detects aflatoxin in the sample. During incubation of the test strip, a continuous band pattern (test line and control line) needs to be formed that is used to determine the concentration of aflatoxin. The intensity of the test line depends on the aflatoxin concentration of the sample.

It increases as the aflatoxin concentration increases. A weak test band can also occur with a negative sample (contaminated below the detection limit). It must be possible to detect the control line after the reaction has taken place in order to verify the functioning of the test. The control line weakens as the aflatoxin concentration in the sample increases.

The test strip is evaluated using the RIDA®SMART APP software (Art. No. ZRSAM).

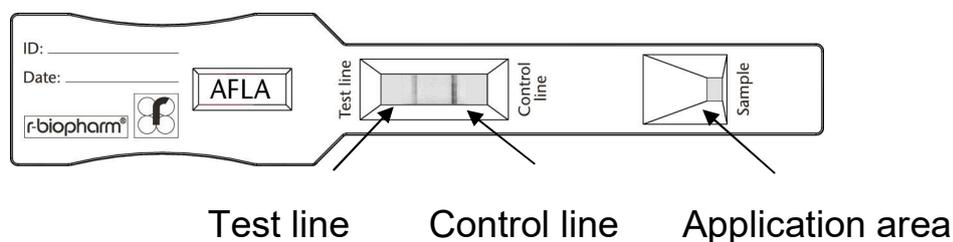
Fig. 1: RIDA®QUICK Aflatoxin RQS test strip.

#### 4. Reagents provided

Each kit contains sufficient materials for max. 20 determinations. Each test kit contains:

Component	Cap color	Format	Volume
20 x Test strip	-	ready to use	One determination per test strip (separately packed)
1 x Extraction solution	green	concentrate 10x	45 mL
1 x Dilution buffer	transparent	ready to use	20 mL
1 x RIDA®SMART APP cover	-	ready to use	Evaluation cover for use with RIDA®SMART APP

#### 5. Equipment and Reagents required but not provided

##### 5.1 Equipment

- Gloves
- Laboratory grinder
- Weighing scale
- Graduated cylinder
- 1 L bottle
- Shaker
- Centrifuge ( $\geq 2000$  g)
- 100  $\mu$ L pipette

- Latest version of the RIDA®SMART APP software (Art. No. ZRSAM) activated and installed on a compatible smartphone.
- Optionally: RIDA®SMART BOX (Art. No. ZRSA-SB)

## 5.2 Reagents

- Distilled (dest.) or deionized (deion.) water

## 6. Warnings and precautions for the users

The product / test is only suitable within the scope of its intended use.

This test should be carried out only by trained laboratory employees. The instruction for use must be strictly followed.

This kit may contain hazardous substances. Please refer to the component safety information in the material safety data sheets (SDS) for this product, available online at [www.r-biopharm.com](http://www.r-biopharm.com).

All reagents and materials must be recovered or disposed after use at customers own responsibility according to the protection of human health and the environment. Please observe the applicable national regulations concerning waste disposal (e.g. Waste Management Act, Regulations on Dangerous Chemicals, etc.).

## 7. Storage instructions

Store the kit at 2 - 8 °C (36 - 47 °F). Do not freeze the test strips or any test kit components.

The test strips are sensitive to moisture. Moist test strips can have a negative impact on the test result. For this reason, it is important to protect the strips against moisture. Do not remove the test strips from the test strip packaging until immediately prior to use in the test.

Do not use the test kit after the expiration date (see test kit label).

Do not interchange individual reagents between kits of different lot numbers.

## 8. Sample preparation

For sample preparation a **ready to use extraction solution** is needed. For this dilute the extraction solution (concentrate 10x) 1:10 (1+9) with distilled water. The extraction solution should be handled under clean working conditions. The ready to

use extraction solution should be stored at 2 - 8 °C (36 - 46 °F) and expires after 8 - 10 weeks. If not all test strips are used within 8 - 10 weeks it is recommended to prepare small volumes of the ready-to-use extraction buffer.

Bring all components necessary for sample preparation to room temperature (20 - 25 °C / 68 - 77 °F) before use in the test and perform sample preparation at room temperature.

The sample should be stored in a cool place, protected against light.

A representative sample (according to accepted sampling techniques) should be ground and thoroughly mixed prior to proceeding with the extraction procedure.

### 8.1 Measurement range 2 - 75 µg/kg

- Weigh 10 g of ground and homogenized sample and add 20 mL of ready-to-use extraction solution
- Shake for 5 min (on shaker)
- Centrifuge (3 min at 2000 g) until a particle-free solution can be obtained
- Use 100 µL of the particle-free solution in the test

#### **Note:**

If necessary, the sample weight can be increased, e.g. 20 g sample in 40 mL ready-to-use extraction solution. In this case the volume of ready-to-use extraction solution has to be adapted.

### 8.2 Measurement range 50 - 300 µg/kg

For an extended measurement range of 50 - 300 µg/kg, further dilute the particle-free solution obtained from sample preparation 8.1.

- Dilute 100 µL of the particle-free solution with 700 µL dilution buffer and mix carefully (= Dilution I)
- Use 100 µL of Dilution I in the test

## **9. Test procedure**

Bring all components necessary for sample preparation to room temperature (20 - 25 °C / 68 - 77 °F) before use in the test and perform sample preparation at room temperature.

- Remove a test strip from its foil pouch.
- Apply 100 µL of prepared sample solution (particle-free solution or dilution I, see 8. Sample preparation) to the application area of the test strip (see Fig. 1).
- Incubate the test strip at room temperature.
- After exactly 5 min of incubation, evaluate the test strip with the RIDA®SMART APP.

## 10. Evaluation

The left line in the reaction area is the test band (test line, see Fig. 1). Its appearance and intensity depends on the aflatoxin concentration of the sample. The sample is contaminated with aflatoxin, if the control band (control line) is visible and the test band (test line) is visible.

The right line in the reaction area is a control band (control line) and must be present after each test procedure. If this line is not visible, the test was not performed correctly or the reagents were not in order. In this case, the test should be repeated with a new test strip. If the right band is missing again, please inform your local distributor.


### Please note:

- Fragmented bands:

Results will be invalid and test strips are inoperative if (testing or control) bands appear fragmented.


- Band detection:

The result field in the RIDA®SMART APP must always be visually checked. The correct test- and control-bands (red bands) have to be framed as shown in Fig. 2. Results are invalid if another band is framed (see e.g. fig. 3).

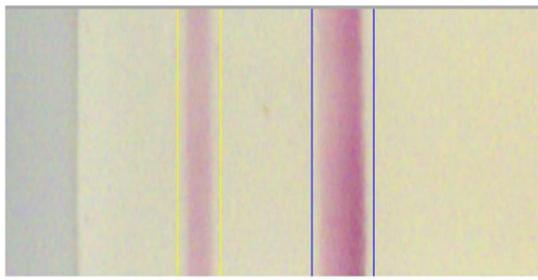Test line Control line

Fig. 2: Band detection valid

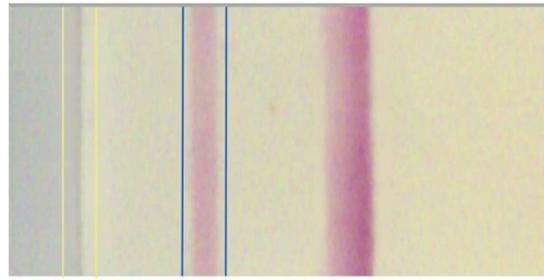Test line Control line

Fig. 3: Band detection invalid

- Reported result:

If the sample was prepared according to 8.1 and is highly contaminated with aflatoxin, the control line will be only faintly visible. In this case the test result might be reported as „>75 ppb or invalid result“ (RIDA®SMART APP). In this instance, dilute the sample as described in 8.2 and repeat the test with a new test strip.

To obtain comparable results, the test strips must always be evaluated after 5 min.

### 10.1 Evaluation with RIDA®SMART APP software in combination with an approved smartphone

First, please read the instructions for use of the RIDA®SMART APP software attentively. The test strip evaluation has to be performed by using the RIDA®SMART APP software according to these guidelines. Please find a short description (Quick Starter Guide) of how to use the RIDA®SMART APP on our website [www.r-biopharm.com](http://www.r-biopharm.com). To analyze this test with RIDA®SMART APP software you must update to the latest software version. On the Certificate of analysis you find the latest version by scanning the QR Code.

For the test strip evaluation the RIDA®SMART APP cover is needed which is included in the test kit.

RIDA®SMART APP software applications for evaluation:

Matrix	Measuring range	RIDA®SMART APP application
Corn	2 – 75 µg/kg	2 – 75 ppb
	50 – 300 µg/kg	50 – 300 ppb

## 10.2 Evaluation with RIDA®SMART APP software in combination with the RIDA®SMART BOX

First, please read the instructions for use of the RIDA®SMART APP software and of the RIDA®SMART BOX attentively. The test strip evaluation with RIDA®SMART APP software in combination with the RIDA®SMART BOX has to be performed according to these guidelines. Please find a short description (Quick Starter Guide) of how to use the RIDA®SMART APP on our website [www.r-biopharm.com](http://www.r-biopharm.com). To analyze this test with RIDA®SMART APP software you must update to the latest software version. On the Certificate of analysis you find the latest version by scanning the QR Code.

The RIDA®SMART software applications that need to be selected can be found in 10.1.

## 11. Recommendation

In order to ensure a high analytical performance we recommend:

- In case of extremely acidic or basic samples, adjust the sample's pH value (pH 6.5 – 7.5) to neutral prior to extraction.
- To do spike experiments to ensure an accurate and correct test procedure.

## 12. Sensitivity

This test determines aflatoxin contaminations in the measurement range of 2 – 300 µg/kg in combination with the RIDA®SMART APP software.

## 13. Further application notes

Further application notes are available on request.

**For further product information and applications, please contact your local distributor or R-Biopharm at this address: [sales@r-biopharm.de](mailto:sales@r-biopharm.de).**

## Version overview

Version number	Chapter and title
2019-11-25	Release version
2020-12-10	General revision
2024-02-19	General revision
2025-02-04	General revision
2025-06-16	Current version Changes made: Art. No. RIDA®SMART APP changed

## Explanation of symbols

- General symbols:

Follow the instructions for use

Batch number

Expiry date (YYYY-MM-DD)

Storage temperature

Article number

Number of test determinations

Manufacturing date (YYYY-MM-DD)

Manufacturer + address

## Disclaimer

1. In conformance with the German Civil Code (“BGB”) R-Biopharm AG provides a limited warranty (“Gewährleistung”) against defects in design and manufacture present at delivery that make the Product unsuitable or unsafe for the contractually specified Product use and location through properly qualified personnel. Such limited warranty warrants Product title and against such material Product defects for 12 months, or if the Product is provided with a shelf life, the length of the declared shelf life, calculated from the date of risk transfer pursuant to delivery terms, and further provided customer provides timely and proper notice of the defect.  
ALL OTHER WARRANTIES OR GUARANTIES, EXPRESS OR IMPLIED, OF ANY KIND ARE EXCLUDED, WHETHER IMPLIED BY CUSTOM, PRACTICE, THE COURSE OF DEALING BETWEEN THE PARTIES OR OTHERWISE. The responsibility for any express or implied warranties or extensions of R-Biopharm’s own warranty made by dealers, distributors, or resellers is solely with such dealer, distributor, or reseller, and R-Biopharm AG expressly disclaims liability for such warranties.
2. R-Biopharm Products are designed for sole use by properly trained and qualified personnel applying appropriate industry standards and practices. Defects are covered only to the extent that such defects are demonstrably the result of defective material, design, parts, or faulty workmanship, which affect safety or result in substandard performance below specifications. R-Biopharm AG is not liable for any improper use nor the consequences of
  - a. the failure to read, understand and follow use and safety instructions, or use proper preparation and storage;
  - b. the failure to utilize trained and unqualified personnel, suitable samples or sampling techniques;
  - c. chemical, electromagnetic, mechanical, or electrolytic influences outside R-Biopharm AG provided standard perimeters through product specific data sheets, written product descriptions, or as otherwise expressly agreed upon in writing; or
  - d. any combination thereof.
3. R-Biopharm AG is also not liable for any changes or modifications, not carried out by R-Biopharm AG, nor consequences of use, applications, or processing which are unsafe or otherwise inconsistent with intended Product purpose as limited by written Product descriptions and specifications of R-Biopharm AG.
4. R-Biopharm AG’s liability for ordinary breach of contract is limited to repair, replacement, other substitute performance, or refund. The choice of remedies is within R-Biopharm AG’s sole discretion. R-Biopharm AG is not contractually liable for any incidental, or consequential damages, including but not limited to purchaser’s expenses, losses, or damages from loss of good will, frustrated business purposes, sales expectations, investment loss, actual or anticipated lost profits, End User indemnity, or any other business expenditures.
5. The foregoing limited warranty is solely intended to fulfill the warranty requirements (“Gewährleistung”) implied by German Civil Code. It is not intended to extend the scope or period of such Gewährleistung or provide additional warranties. Nor is this Disclaimer otherwise intended to limit or extent mandatory liability for damages in tort resulting from injury to life, body, or health, for intentional or grossly negligent acts, fraud, or due to strict liability under applicable product liability laws. A reversal of the burden of proof or rules of contract interpretation is not intended.

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are registered trademarks of R-Biopharm AG.