



R-BIOPHARM
RHÔNE LTD

Accessory products – for your convenience



Complements immunoaffinity product range




Technical support available from our team



Technical videos available on YouTube

More information:

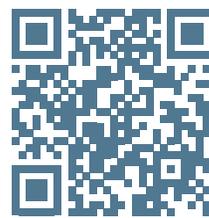

<https://r-b.io/food>

Accessory products – for your convenience

Immunoaffinity column Accessory Pack

Art. no.: AP01

- Contains 10 glass barrels, 10 syringe pump units and 10 reusable PTFE adapters for use with immunoaffinity columns.
- An empty glass barrel can be stacked on top of an immunoaffinity column to act as a larger load reservoir.
- Combine two immunoaffinity columns to provide different selectivities.
- The PTFE adapters can be used in positive pressure methods.
- Can be used in conjunction with immunoaffinity columns to aid the analysis.
- All components can be re-used following cleaning.
- Also suitable for use with OCHRACARD.

Immunoaffinity column rack

Art. no.: CR1 / CR2

- The simple and easy-to-use racking system has been developed for use when performing immunoaffinity column analysis. Inserts are also available for the rack.
- Increased sample productivity due to multiple column holders.

- Lightweight structure allow easy movement around the laboratory.
- Chemical resistant. The rack will keep its appearance for longer.
- Rack dimensions:
 - Overall (45 cm long x 15 cm deep x 40 cm high).
 - Holes (6, each 18 mm diameter)
- Easy to clean and decontaminate.
- Height adjustable to a maximum of 40 cm.
- Inserts (CR2) are available allowing it to be used with vials or evaporation tubes.
- Low cost.

PBS Tablets

Art. no.: RP202

Phosphate-buffered saline (PBS) is a buffer solution used in biological research. It is a water-based salt solution containing sodium phosphate, sodium chloride, potassium chloride and potassium phosphate. The osmolality and ion concentrations of the solutions match those of the human body (isotonic) and are non-toxic to most cells.

Ordering information

Product	Description	No. of tests / amount	Art. No.
	for immunoaffinity columns		
PBS tablets	Phosphate Buffered Saline Tablets	100 (suitable for 10 L)	RP202
Immunoaffinity column rack	Durable brass and PFTE rack allowing 6 samples to be processed at one time using immunoaffinity columns	1 unit	CR1
Immunoaffinity columns accessory pack	Glass barrels, syringes and adapters for use with all formats of R-Biopharm Rhône immunoaffinity columns	10 each	AP01