



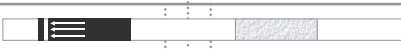


bioavid Hook-Linie – Auswertekarte

bioavid Lateral Flow Testkits mit Hook-Linie – BLH-Serie


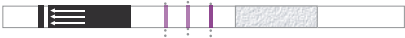


Inkubationszeit: 5 min

	T-Linie Testlinie	H-Linie Hook-Linie	C-Linie Kontrolllinie	
Negativ	-	+	+	
Positiv	+	+	+	
Positiv	+	+/-	+	
Hoch Positiv höhere Verdün- nung testen	-	-	+	
Invalide	+/-	+/-	-	

bioavid hook line – evaluation card

bioavid lateral flow test kits with hook line – BLH-series

Incubation time: 5 min

	T-line test line	H-line hook line	C-line control line	
Negative	-	+	+	
Positive	+	+	+	
Positive	+	+/-	+	
High Positive higher dilution to be tested	-	-	+	
Invalid	+/-	+/-	-	